P/ACE™ MDQ Plus Capillary Electrophoresis System

System Overview Guide

RUO-IDV-05-4901-A
January 2019

This document is provided to customers who have purchased SCIEX equipment to use in the operation of such SCIEX equipment. This document is copyright protected and any reproduction of this document or any part of this document is strictly prohibited, except as SCIEX may authorize in writing.

Software that may be described in this document is furnished under a license agreement. It is against the law to copy, modify, or distribute the software on any medium, except as specifically allowed in the license agreement. Furthermore, the license agreement may prohibit the software from being disassembled, reverse engineered, or decompiled for any purpose. Warranties are as stated therein.

Portions of this document may make reference to other manufacturers and/or their products, which may contain parts whose names are registered as trademarks and/or function as trademarks of their respective owners. Any such use is intended only to designate those manufacturers' products as supplied by SCIEX for incorporation into its equipment and does not imply any right and/or license to use or permit others to use such manufacturers' and/or their product names as trademarks.

SCIEX warranties are limited to those express warranties provided at the time of sale or license of its products and are SCIEX's sole and exclusive representations, warranties, and obligations.

SCIEX makes no other warranty of any kind whatsoever, expressed or implied, including without limitation, warranties of merchantability or fitness for a particular purpose, whether arising from a statute or otherwise in law or from a course of dealing or usage of trade, all of which are expressly disclaimed, and assumes no responsibility or contingent liability, including indirect or consequential damages, for any use by the purchaser or for any adverse circumstances arising therefrom.


For Research Use Only. Not for use in Diagnostic Procedures.

AB Sciex is operating as SCIEX.

The trademarks mentioned herein are the property of AB Sciex Pte. Ltd. or their respective owners.

AB SCIEX™ is being used under license.

© 2019 AB Sciex

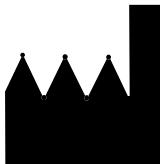AB Sciex LLC
500 Old Connecticut Path
Framingham, Massachusetts 01701
USA

Contents

Chapter 1 Safety, Notices and Labels	7
Overview	7
Alerts for Warning, Caution, Important, and Note	8
Instrument Safety Precautions	9
Moving Parts or Sharp Objects	10
Electrical Safety	10
Laser Safety (for Optional Laser Device)	10
Class 1 Laser Caution Label	11
Chemical Precautions	11
Safety Symbols and Labels	12
High Voltage Electric Shock Risk Symbol	12
Attention Safety Symbol	12
Sharp Object Label	12
Fuse Warning Symbol	13
Cancer and Reproductive Harm Label	13
RoHS Notices	13
China RoHS Caution Label	13
Other Instrument Labels	14
Recycling Label (WEEE)	14
Disposal of Devices Containing Mercury Components	14
CE Mark Label	15
RCM Mark Label	15
CSA Mark Label	15
Chapter 2 安全、通知和标签	17
综述	17
警告、注意、重要事项以及注释的提示	18
仪器安全防护措施	19
活动部件或尖锐物	20
电气安全	20
激光安全 (适用于可选激光设备)	20
1类激光小心标签 (1级激光警告标签)	20
化学品注意事项	21
安全符号和标签	21
高压电击危险符号	21
注意安全符号	21
尖锐物标签	22
保险丝警告符号	22
致癌和生殖危害标签	22

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How to Load the Cartridge and Samples	52
How to Run a Method.....	54
Analysis	57
Area Calculation	60
Annotations	61
Peak Identity	62
Run Multiple Methods with a Sequence.....	62
Viewing Reports	66
Standard Reports	66
Method Custom Reports	66
Sequence Custom Report	66
Advanced Reports	66
Suitability Reports.....	67
Viewing and Printing Reports.....	68
Maintenance	68
Appendix A Specifications	69
P/ACE™ MDQ Plus System	69
Validated Controller Configuration	70
Sample Temperature Control	70
Capillary Temperature Control.....	70
Pressure and Vacuum System.....	71
Detector Specifications.....	71
UV Detector	71
Laser Induced Fluorescence (LIF) Detector (Optional)	72
Photo Diode Array (PDA) Detector (Optional)	72
Contact Us	75
Customer Training	75
Online Learning Center	75
SCIEX Support.....	75
CyberSecurity	75
Documentation	76

Safety, Notices and Labels

Overview

A description of the symbols and labels that are used on the SCIEX P/ACE™ MDQ Plus Capillary Electrophoresis System, or that are shown in this manual, can be found in this section.

Do not attempt to perform any procedure before carefully reading all instructions. If in doubt as to how to proceed in any situation, contact your SCIEX representative.

SCIEX urges its customers and employees to comply with all national health and safety standards such as the use of barrier protection. This may include, but is not limited to, protective eyewear, gloves, and suitable laboratory attire when operating or maintaining this or any other automated laboratory instrumentation.

If the equipment is used in a manner not specified by SCIEX, the protection provided by the equipment may be impaired.

Alerts for Warning, Caution, Important, and Note

All Warnings and Cautions in this document include an exclamation point, framed within a triangle.

The exclamation point symbol is an international symbol which serves as a reminder that all safety instructions should be read and understood before installation, use, maintenance, and servicing are attempted.

WARNING indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury. It may also be used to indicate the possibility of erroneous data that could result in an incorrect diagnosis. May be used to indicate the possibility of severe instrument damage.

CAUTION indicates a potentially hazardous situation, which, if not avoided, may result in minor or moderate injury. It may also be used to alert against unsafe practices. May be used to indicate the possibility of erroneous data that could result in an incorrect diagnosis.

IMPORTANT Used for comments that add value to the step or procedure being performed. Following the advice in the Important adds benefit to the performance of a piece of equipment or to a process.

NOTE Used to call attention to notable information that should be followed during installation, use, or servicing of this equipment.

Instrument Safety Precautions

WARNING

Risk of operator injury if:

- All doors, covers and panels are not closed and secured in place prior to and during instrument operation.
- The integrity of safety interlocks and sensors is compromised.
- You contact moving parts.
- You mishandle broken parts.
- Doors, covers and panels are not opened, closed, removed and/or replaced with care.
- Improper tools are used for troubleshooting.

To avoid injury:

- Keep doors, covers and panels closed and secured in place while the instrument in use.
- Take full advantage of the safety features of the instrument. Do not defeat safety interlocks and sensors.
- Acknowledge and act upon instrument alarms and error messages.
- Keep away from moving parts.
- Report any broken parts to your SCIEX representative.
- Use the proper tools when troubleshooting.

CAUTION

System integrity could be compromised and operational failures could occur if:

- This equipment is used in a manner other than specified. Operate the instrument as instructed in the product manuals.
- You introduce software that is not authorized by SCIEX into your computer. Only operate your system's computer with software authorized by SCIEX.
- You install software that is not an original copyrighted version. Only use software that is an original copyrighted version to prevent virus contamination.

CAUTION

If you purchased this product from anyone other than SCIEX or an authorized SCIEX distributor, and, if it is not presently under a SCIEX Service Maintenance Agreement, SCIEX cannot guarantee that the product is fitted with the most current mandatory engineering revisions or that you will receive the most current information bulletins concerning the product. If you purchased this product from a third party and would like further information concerning this topic, contact your

SCIEX representative.

Moving Parts or Sharp Objects

Risk of personal injury. To avoid injury due to moving parts, observe the following:

- Never attempt to exchange labware, reagents, or tools while the instrument is operating.
- Never attempt to physically restrict any of the moving components of the instrument.
- Keep the instrument work area clear to prevent obstruction of the movement.

Electrical Safety

To prevent electrically related injuries and property damage, properly inspect all electrical equipment prior to use and immediately report any electrical deficiencies. Contact a SCIEX representative for any servicing of equipment requiring the removal of covers or panels.

Laser Safety (for Optional Laser Device)

This product may contain a laser module. The laser (optional) is designated as "Class 3B." The "3B" classification means that "direct intrabeam viewing of this type of laser is always hazardous to personnel."

The laser and several other integral components are contained in a sealed housing that together comprise the laser assembly. The laser assembly has no user serviceable parts. Service of the laser assembly is restricted to qualified SCIEX Field Service Employees (FSE).

Therefore, the overall laser classification of the CE Instrument is "Class 1," defined as "lasers which are safe under reasonably foreseeable conditions of operation."

To prevent users from potentially harmful laser light, observe all safety warnings and NEVER REMOVE THE OUTER CASING OF THE LASER ASSEMBLY.

IMPORTANT The laser markings noted above can vary depending on the type of device, and will be described in the documentation provided with the module.

Class 1 Laser Caution Label

If the instrument contains a laser system, a label reading “THIS PRODUCT CONFORMS TO APPLICABLE REQUIREMENTS OF 21 CFR 1040 AT THE DATE OF MANUFACTURE” is found near the Name Rating tag. The laser light beam is not visible.

Chemical Precautions

- Determine which chemicals have been used in the system prior to service and regular maintenance. Refer to Safety Data Sheets for the health and safety precautions that must be followed with chemicals.
- Work in a well-ventilated area.
- Always wear assigned personal protective equipment, including powder-free neoprene or nitrile gloves, safety glasses, and a laboratory coat.
- Follow required electrical safe work practices.
- Avoid ignition sources when working with flammable materials, such as isopropanol, methanol, and other flammable solvents.
- Take care in the use and disposal of any chemicals. Potential risk of personal injury if proper procedures for handling and disposing of chemicals are not followed.
- Avoid skin contact with chemicals during cleaning and wash hands after use.
- Comply with all of the local regulations for the storage, handling, and disposal of biohazardous, toxic, or radioactive materials.
- The lamp in this product contains mercury. Do not put in the trash. Recycle or dispose of according to local, state, or federal laws.

Safety Symbols and Labels

High Voltage Electric Shock Risk Symbol

This symbol indicates that there is high voltage and there is a risk of electric shock, and the operator should use care when accessing this area.

Attention Safety Symbol

This symbol calls attention to important information to read, or is accompanied by another symbol indicating a particular safety hazard. The information is located either on the label with the symbol or in the P/ACE™ MDQ Plus documentation.


Sharp Object Label

This symbol indicates that there are sharp objects, and the operator should use care when accessing this area.

Fuse Warning Symbol

This symbol indicates that the instrument fuse may be replaced only with the specified type and rating.

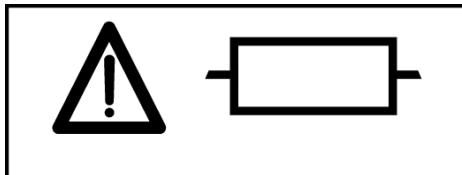Cancer and Reproductive Harm Label

This label indicates there is risk of cancer and reproductive harm to the operator. This warning is called for by the California Safe Drinking Water and Toxic Enforcement Act of 1986, commonly referred to as Proposition 65, enacted by the State of California. A full list of harmful chemicals is available at www.P65Warnings.ca.gov.

RoHS Notices

These labels and materials declaration table (the Table of Hazardous Substance's Name and Concentration) are to meet People's Republic of China Electronic Industry Standard SJ/T11364-2006 “Marking for Control of Pollution Caused by Electronic Information Products” requirements.

China RoHS Caution Label

This label indicates that the electronic information product contains certain toxic or hazardous substances. The center number is the Environmentally Friendly Use Period (EFUP) date, and indicates the number of calendar years the product can be in operation. Upon the expiration of the EFUP, the product must be immediately recycled. The circling arrows indicate the product is recyclable. The date code on the label or product indicates the date of manufacture.

Other Instrument Labels

This section provides information for some labels and symbols appearing on the P/ACE™ MDQ Plus instrument housing. These labels and symbols may be associated with user-serviceable procedures. Individual hazards associated with a specific procedure in this manual may use these labels and symbols, and are included in Warnings or Cautions within the procedures for that task.

Recycling Label (WEEE)

The symbol of a crossed-out wheeled bin on the product is required in accordance with the Waste Electrical and Electronic Equipment (WEEE) Directive of the European Union.

The presence of this marking on the product indicates:

- That the device was put on the European Market after August 13, 2005 and
- That the device is not to be disposed via the municipal waste collection system of any member state of the European Union.

Follow local municipal waste ordinances for proper disposal provisions to reduce the environmental impact of WEEE (waste, electrical, and electronic equipment). To safely dispose of this equipment, contact a local Customer Service office for complimentary equipment pick-up and recycling.

Disposal of Devices Containing Mercury Components

This product may contain a mercury-added part. Recycle or dispose of according to local, state, or federal laws. It is very important that you understand and comply with the safe and proper disposal of devices containing mercury components (switch, lamp, battery, relay, or electrode). The mercury component indicator label can vary depending on the type of device.

The lamp in this product contains mercury. Do not put in the trash. Recycle or dispose of according to local, state or federal laws.

CE Mark Label

A “CE” mark indicates that a product has been assessed before being placed on the market, and has been found to meet European Union safety, health, and/or environmental protection requirements.

RCM Mark Label

The RCM mark is intended for use on products that comply with Australian communications Media Authority (ACMA) EMC Requirements.

CSA Mark Label

The CSA symbol on the P/ACE™ MDQ Plus instrument indicates that the instrument has been certified by a Nationally Recognized Testing laboratory (NRTL) to applicable Laboratory Equipment Safety Standards for the United States and Canada.

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
System Overview

Instrument


The main components of the SCIEX P/ACE™ MDQ Plus Capillary Electrophoresis System include trays that hold vials of sample, buffer, and other solutions, an interface block, a high-voltage power supply and electrodes, a source optics module and detector, temperature control hardware, and a sample injection mechanism (refer to [Figure 4.1](#), [Figure 4.2](#), and [Figure 4.3](#)).

Figure 4.1 P/ACE™ MDQ Plus Instrument—Front View

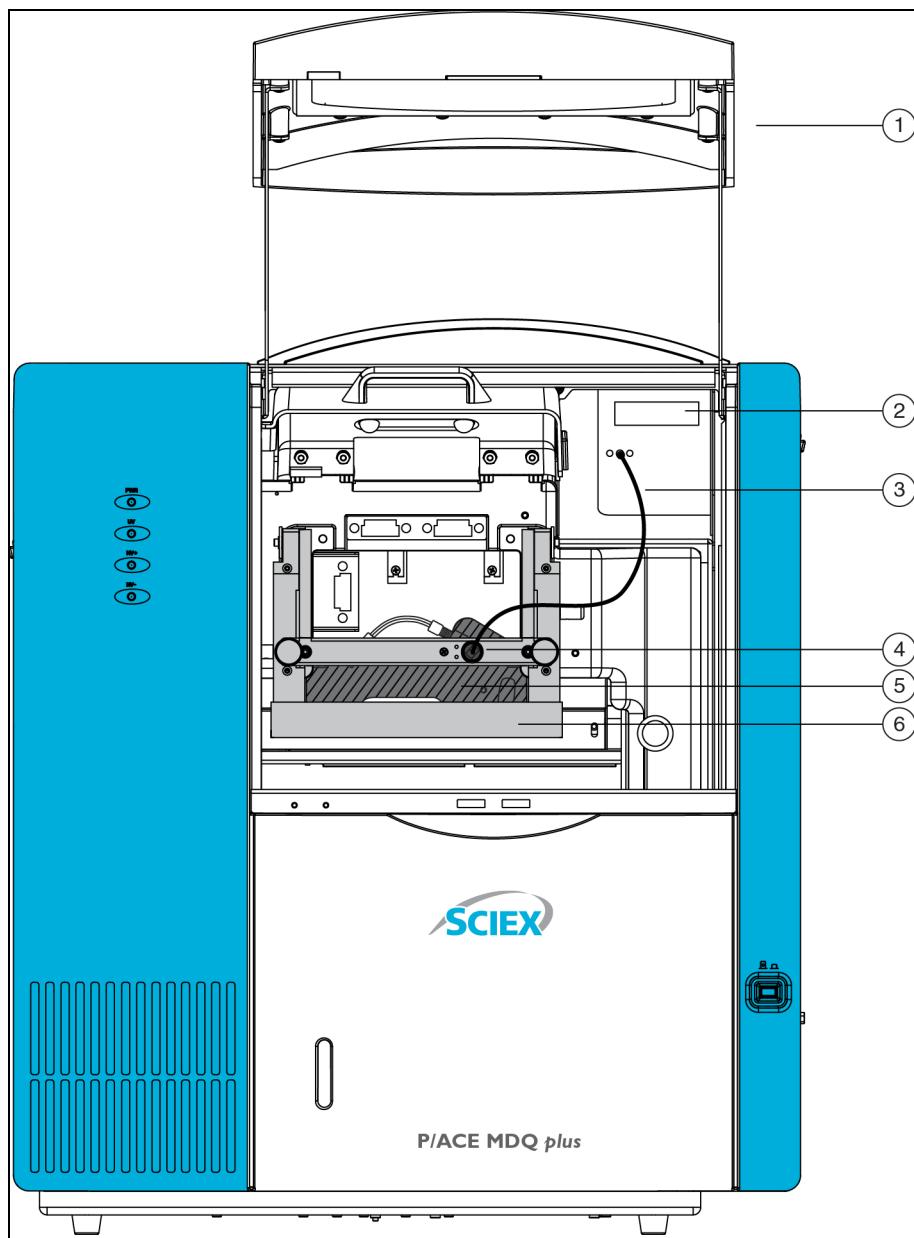- | | |
|--------------------------------|-----------------------------------|
| 1. Outer Door (open position) | 4. Clamp Bar and Cable Connection |
| 2. UV Detector | 5. Capillary Cartridge |
| 3. Two-ended Fiber Optic Cable | 6. Interface Block |

Figure 4.2 P/ACE™ MDQ Plus Instrument—Side View

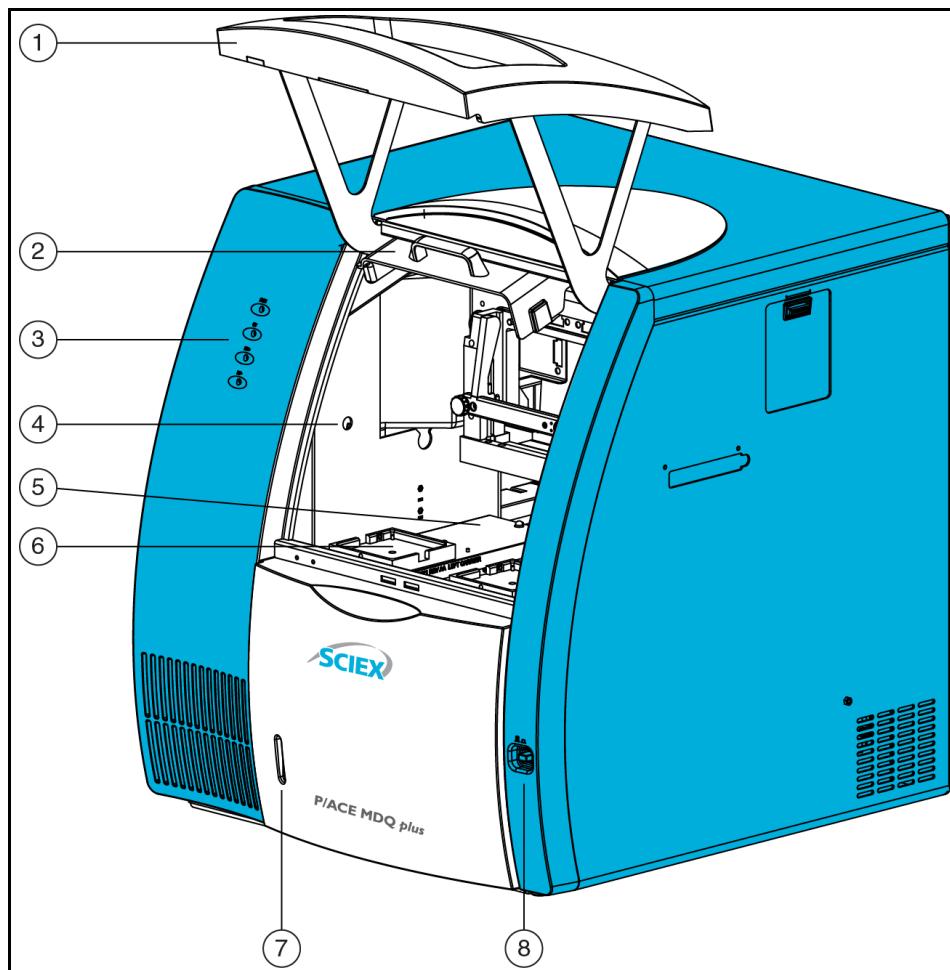
- | | |
|---|---------------------------|
| 1. Outer Door or Sample Cover (open) | 5. Sample Trays |
| 2. Inner Door or Cartridge Cover (open) | 6. Buffer Trays |
| 3. Indicators | 7. Fluid Bubble Indicator |
| 4. Fluid Fill Port | 8. Convenience Switch |

The convenience switch is on the lower-right side of the front of the instrument. All connections for external system components are on the upper-left side panel of the instrument, except for the AC inlet and the fuse holder. Three fans supply cooling air flow for internal system components. Air is exhausted through the vents at the side and back of the instrument. Keep at least six inches of clearance at each vent to ensure adequate air flow.

Sample Handling System

The sample handling system holds four trays; two sample trays (inlet and outlet), and two buffer trays (inlet and outlet). The sample trays are primarily used for samples; the buffer trays hold the buffer and rinse solutions required for electrophoresis. The trays are on two parallel tracks. Under normal operating conditions, the trays on the left are inlet trays for sample and buffer; the trays on the right are outlet trays for sample and buffer ([Figure 4.3](#)).


Figure 4.3 P/ACE™ MDQ Plus Trays

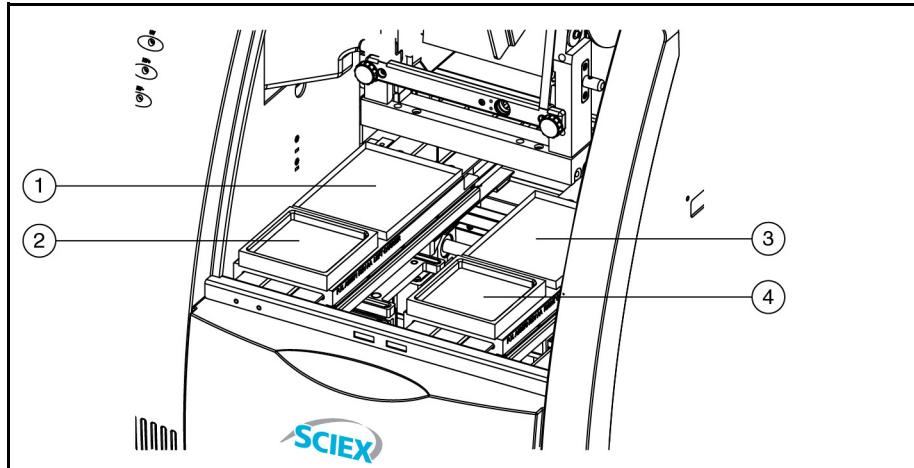- | | |
|--|---|
| 1. Inlet Sample Tray (48) or 96-Well Plate | 3. Outlet Sample Tray (48) or 96-Well Plate |
| 2. Inlet Buffer Tray (36 Vials) | 4. Outlet Buffer Tray (36 Vials) |

Each buffer tray has slots for 36 universal vials. Sample trays hold either 48 universal vials or a 96-well plate. Each slot is assigned a number from the front to the back, starting with the number 1, and assigned a letter from left to right, starting with the letter A ([Figure 4.4](#)).

Figure 4.4 Sample and Buffer Trays

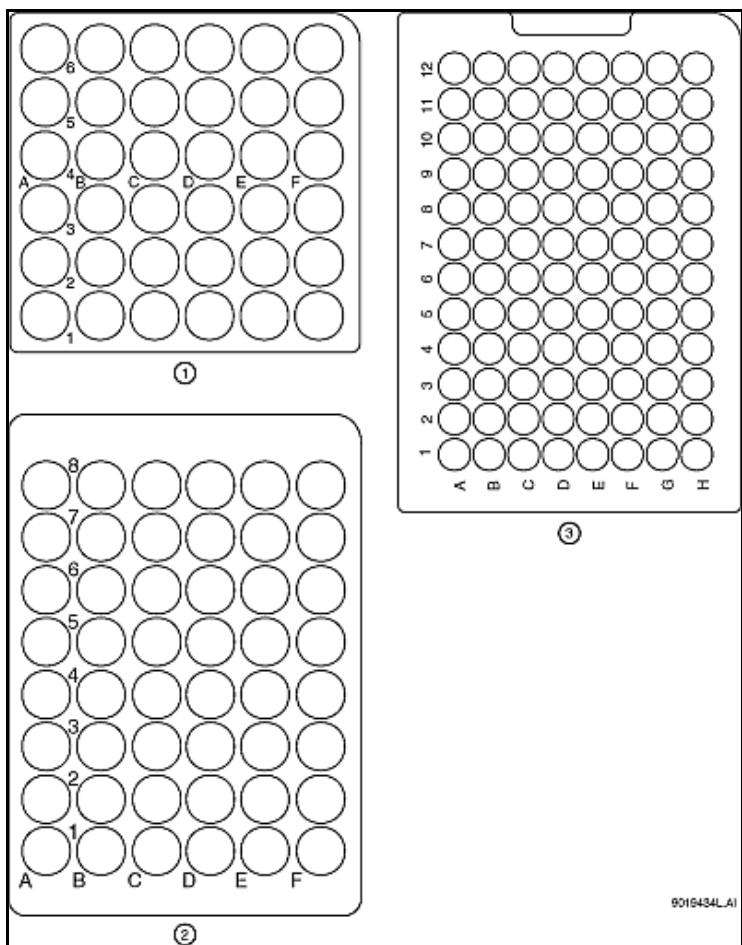1. Buffer Tray
2. 48-Vial Sample Tray for universal vials. Universal vials can also be used as a holder for micro vials containing samples.
3. 96-Position Sample Tray

⚠️ WARNING

The P/ACE™ MDQ Plus instrument is not intended to accommodate volatile materials in 96-well plates. Volatile solvents can release hazardous or flammable vapors leading to a risk of fire or explosion. The solvent vapors can damage the instrument. Do not use volatile solvents in 96-well plates.

⚠️ WARNING

Wear safety glasses when opening the sample cover while the vials are pressurized.


The universal vials are pressurized during rinse and separation-with-pressure events. To reduce the risk of breakage and expelled particles, use only universal plastic vials (PN A62251), and inspect every vial for

damage before use. Do not use any vial that appears cracked or damaged in any way.

Capillary Cartridge

The separation capillary is installed in a cartridge. The cartridge design protects the capillary, supplies a path for liquid coolant, simplifies installation into the instrument, and aligns the detection window in the optics. The components of the cartridge are shown in [Figure 4.5](#).

Figure 4.5 The Capillary Cartridge

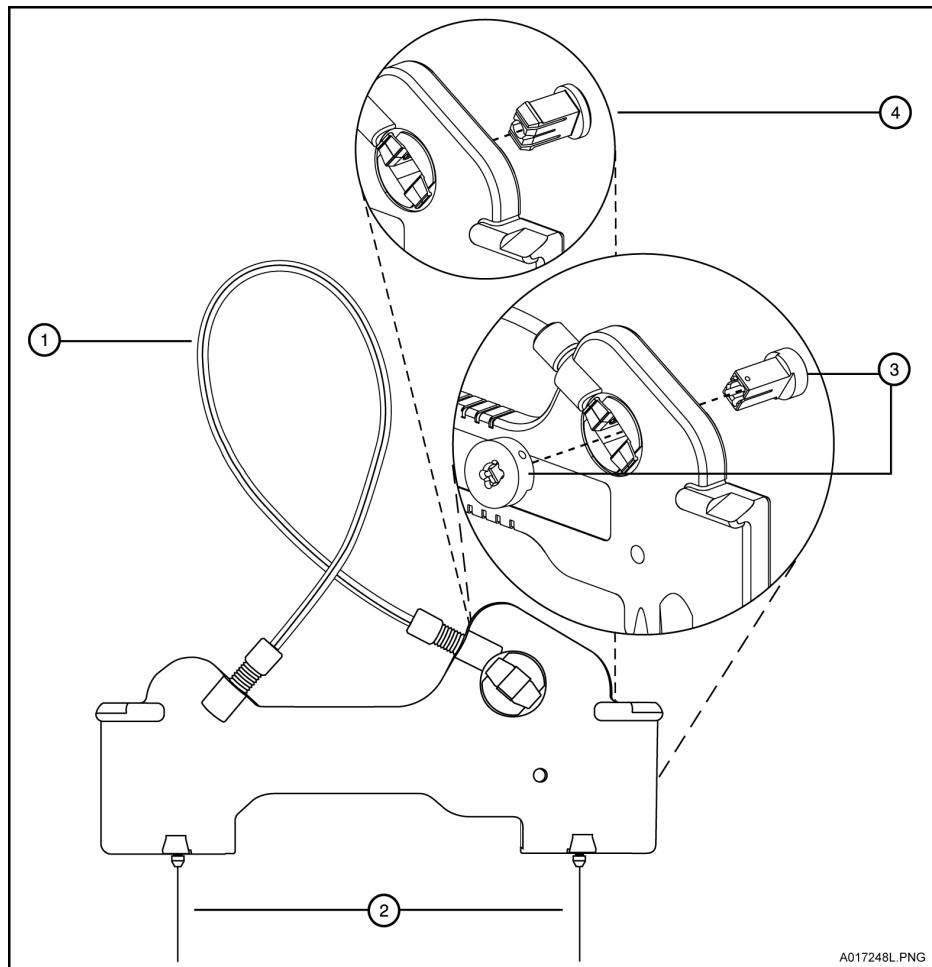1. Coolant tubing with capillary inside
2. Double seal
3. Detector Window and Aperture - LIF Detector
4. Detector Window and Aperture - UV and PDA Detectors

The detection window is an area of the capillary where the polyimide coating is removed to show the transparent fused-silica. This area of the capillary is put in a part of the cartridge that contains a plug that connects the window to the optical system. One type of plug is used for UV and PDA detectors, and a second type is used for LIF detectors. The procedures to install a capillary in a cartridge are included in the *System Maintenance Guide*.

The capillary temperature is controlled with an inert liquid that circulates through the cartridge. The temperature is controlled in a range from 10 °C below ambient (with a minimum of 15 °C) to 60 °C. Coolant flows through the cartridge through two openings in the bottom of the housing (found between the ends of the capillary). This fluid removes the heat generated by electrophoresis.

Fluid Delivery, Power Supply, and Interlock

Syringe Pump

The P/ACE™ MDQ Plus instrument can generate pressures with an internal pump mechanism. This pump can supply 0.1 psi to 25 psi to perform pressure injections or low-pressure mobilizations. The pump can apply a maximum of 100 psi to move the fluids through the capillary. Vacuum injections can be performed from 0.1 psi to 5.0 psi. The pressure can be applied to both ends of the capillary at the same time to prevent outgassing of gels.

High Voltage Power Supply

The high-voltage power supply can deliver a maximum of 30 kV with a maximum current of 300 μA. The voltage range is from 1 kV to 30 kV in 100 V increments. The polarity is configured in the software. The current range is from 3.0 μA to 300 μA in 0.1 μA increments. The software allows the user to select current, voltage, or power operation. During operation, the system will ramp the voltage or current up to the programmed value. Limits for voltage, current and power can be entered to protect the capillary. For example, if the user programs a voltage setting for 30 kV, but the setting for current is only 3.0 μA, the system can reach the limit set for current before reaching the voltage setting, and control voltage to keep that current.

LED Indicators

The front panel of the instrument contains LED indicators for power, UV, and high voltage ([Figure 4.2](#)).

Cartridge and Sample Cover Interlocks

The hinged doors of the P/ACE™ MDQ Plus instrument have interlock sensors that prevent unsafe access to the inside of the instrument. The first door is called the outer door or sample cover; the second is called the inner door or cartridge cover ([Figure 4.2](#)).

Opening the sample cover:

- Stops any tray movement immediately.
- Prevents the execution of any programmed events that require tray movement.
- Aborts a method when a step that requires tray movement is encountered.

Opening the cartridge cover:

- Turns off the high voltage if it is on.
- Turns off the pump that circulates the capillary coolant.
- Moves the detector filter wheel to the closed position.

UV Detector Optics

The UV optics include an ultraviolet light source, wavelength filters, aperture, capillary, and a photodiode detector, as shown in [Figure 4.6](#).


The UV source is a deuterium lamp with a wavelength range of 190 nm to 600 nm. Two lenses focus and direct the output of the lamp through one of the wavelength-selecting filters found in a rotating wheel behind the capillary cartridge. The beam continues through the aperture in the cartridge plug and through a section of the capillary that has been treated to remove the polyimide coating (the detection window). The non-absorbed beam then continues through a fiber optic cable to a photodiode. The light signal is converted to an electrical signal, digitized, and sent to the 32 Karat workstation for processing by the software. This signal is also available as an analog output through a connection on the left side of the instrument.

The design of the instrument insures that the optical system stays in alignment. No user alignments are required.

There are eight positions on the UV filter wheel. UV detector systems are shipped with four standard filters: 200 nm, 214 nm, 230 nm, and 254 nm (10 nm bandwidth). The filters are installed in positions 2, 3, 4 and 5, respectively, on the filter wheel. Position 1 is opaque and functions as a shutter for the detection system.

Additional wavelengths are obtained by placing the appropriate filters in positions 6, 7, and 8. If required the standard filters can be replaced. If the instrument is used with a PDA detector, position 8 must be left open (no filter). The filter wheel will accommodate $\frac{1}{2}$ inch (12.7 mm) diameter filters with wavelengths from 190 nm to 600 nm.

Figure 4.6 UV Optics Layout

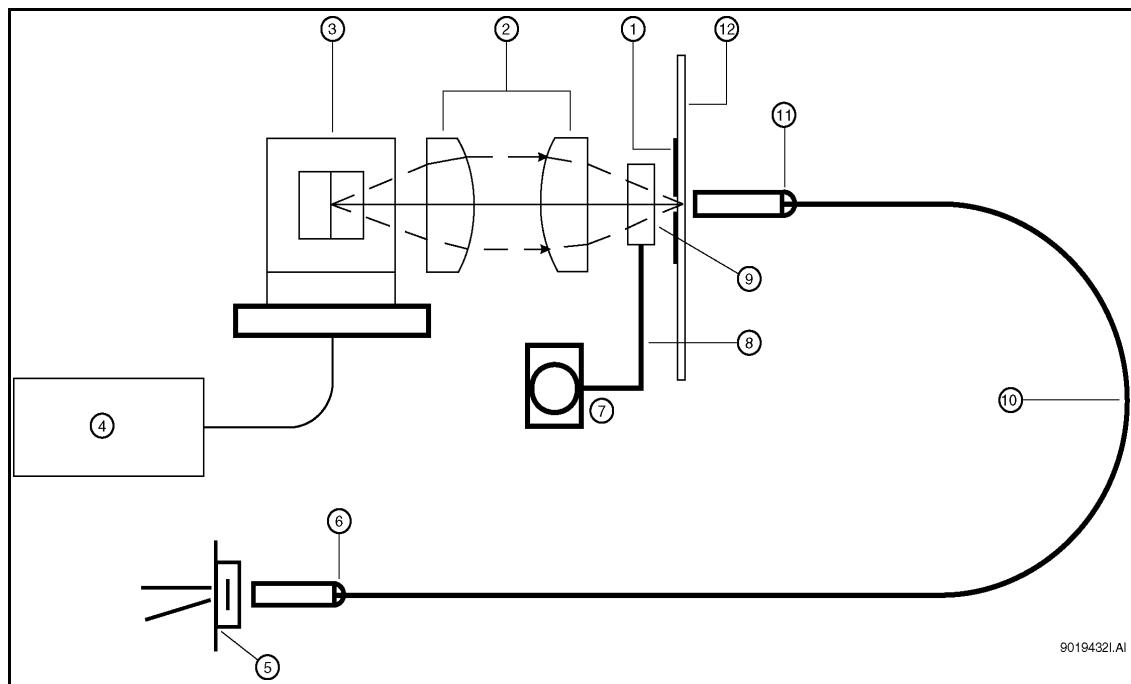- 1. Capillary Aperture
- 2. Lenses
- 3. Deuterium Lamp
- 4. Lamp Power Supply
- 5. Photodiode
- 6. Fiber Optic Connection
- 7. Motor
- 8. Position Filter Wheel
- 9. Filter Position (for example, 214 nm)
- 10. Fiber Optic Cable
- 11. Fiber Optic Connector
- 12. Capillary

The Photo Diode Array (PDA) Detector (Optional)

The optional photo diode array (PDA) detector, like the UV detector, uses the absorbance of light to detect if there are samples as they go through the detection window. Unlike the UV detector, the PDA detector can give spectral analysis of samples.

The PDA detector uses the same cartridge configuration as the UV detector. Refer to the *System Maintenance Guide* for a description of the cartridge.

In PDA detection, the full spectrum of light from the deuterium lamp illuminates the capillary (Figure 4.7). Light that is not absorbed by samples is delivered by a fiber optic cable to a grating that breaks the light into a spectrum. This spectrum is projected onto an array of 256 photodiodes. With this design, the absorbance profile of the sample is measured. The PDA detector also allows the simultaneous measurement of light at different discrete wavelengths. The photo diode array converts the light signal into an electrical signal. The electrical signal is digitized and sent to the 32 Karat workstation for processing by the software.


The PDA detector always uses filter wheel position 8. When the PDA is in use, it is essential that no filter is present in position 8.

The PDA detector is calibrated using discrete emission wavelength bands generated by a mercury lamp. The mercury lamp is an important part of the detector system. When requested by the user, the calibration is performed automatically.

⚠️ WARNING

Do not put the mercury lamp into the regular trash. Mercury is a hazardous material and must be disposed of in accordance with local, state, and federal laws.

Figure 4.7 Diode Array Optics Layout

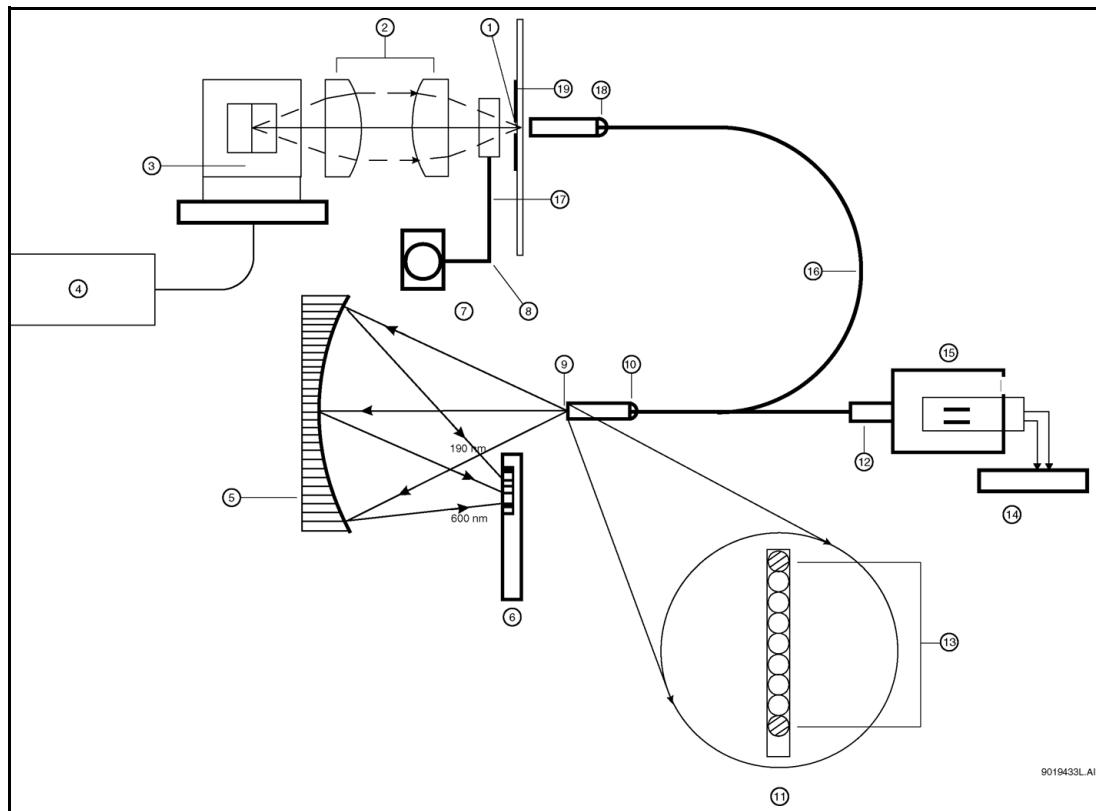- | | |
|--------------------------------|---|
| 1. Capillary Aperture | 11. 9 μm by 200 μm Fiber Array (Slit) |
| 2. Lenses | 12. Fiber Optic Connector |
| 3. Deuterium Lamp | 13. Mercury Calibration Fibers |
| 4. Lamp Power Supply | 14. Mercury Lamp Power Supply |
| 5. Concave Holographic Grating | 15. Mercury Lamp |
| 6. 256 Element Diode Array | 16. Y-fiber Optic Cable |
| 7. Motor | 17. Filter wheel in open position 8 |
| 8. Position Filter Wheel | 18. Fiber Optic Connector |
| 9. Monochromator Entrance Slit | 19. Capillary |
| 10. Fiber Optic Connector | |


The Laser Induced Fluorescence (LIF) Detector (Optional)

The optional laser induced fluorescence (LIF) detector consists of the LIF detector module and a laser module. A capillary cartridge with an LIF detector plug installed is required for use with this system.

The LIF detector uses an integrated 488 nm laser light source. A fiber cable transmits excitation light from the laser to the capillary in the cartridge. Substances in the capillary that fluoresce at the laser wavelength are detected. The LIF detector measures and records this fluorescence, which appears as a peak on the electropherogram. For more information on the LIF system, refer to the *System Maintenance Guide*.

The initial installation of the LIF detector is performed by a SCIEX Field Service Employee (FSE). The P/ACE™ MDQ Plus instrument can easily be changed between UV and the optional LIF and PDA modes because the detector components are modular. Refer to [Figure 4.8](#) for a stylized diagram of the LIF optical system.

Figure 4.8 LIF Optical System

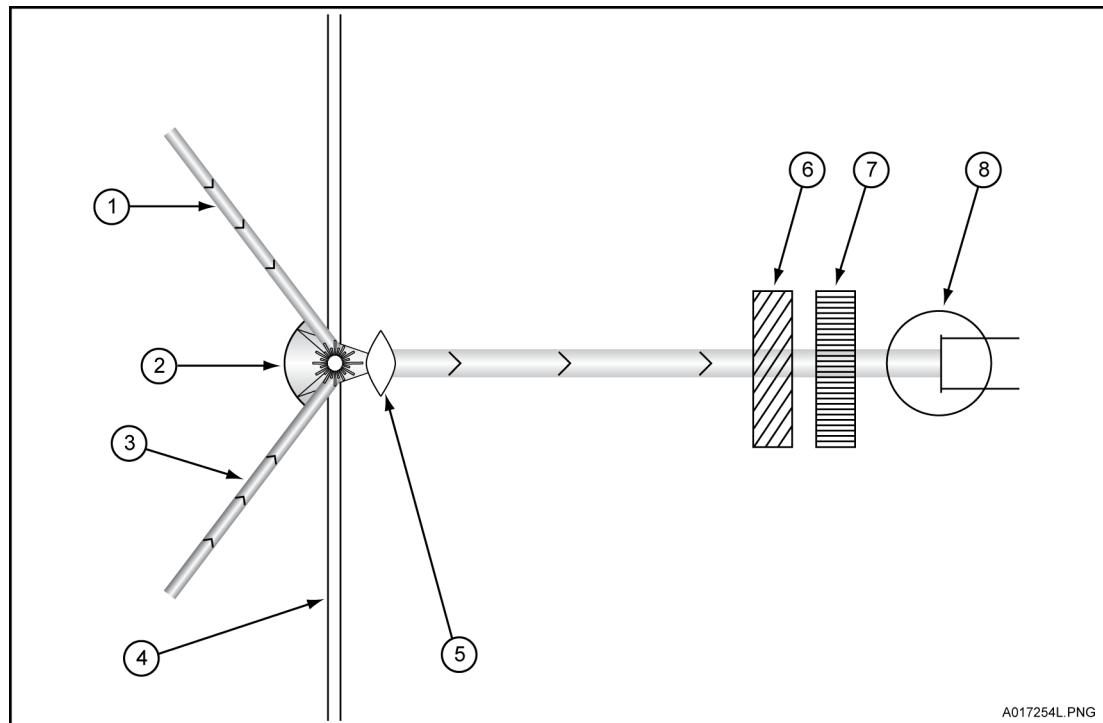- | | |
|---------------------|--------------------------|
| 1. Laser 1 | 5. Ball Lens |
| 2. Spherical Mirror | 6. Laser Filters |
| 3. Laser 2 | 7. Emission Filters |
| 4. Capillary | 8. Photo Multiplier Tube |

Laser Module (Optional)

The following section describes the 488 nm laser module and how it interfaces with the P/ACE™ MDQ Plus instrument and the LIF detector.

 **WARNING**

During normal operation of the LIF detector, laser light is not accessible to the user. To prevent potentially harmful laser light from being emitted from the end of the fiber cable, an interlock mechanism turns off the laser if the laser fiber optic cable is disconnected from the interconnect module or if the cartridge cover is opened.

Always turn off the P/ACE™ MDQ Plus instrument before removing any of the LIF system modules.

488 nm Laser Module

The 488 nm laser is a solid state laser mounted inside the instrument. It has a fiber connector and an electrical interlock connector, for an external laser connection, that is accessible from the right side of the instrument. This external laser connection accommodates a fiber coupled external laser source for LIF detection.

CHAPTER 5

System Procedures and Training Guide

Controller and Instrument Start Up

Controller/Network Logon

The P/ACE™ MDQ Plus software, which includes the 32 Karat software, runs on a Windows 10 controller. If the instrument is installed on a network, the user name and password will be supplied by the network administrator. The 32 Karat software functionality is different when it is run on a network versus as a stand-alone system. In particular, the system administration features can use network names and passwords automatically. The default network identification is **32 Karat** and the work group name is **WORKGROUP**. If more than one P/ACE™ MDQ Plus system is going to be installed on the network, unique network names are required for each workstation. Normally these changes are performed by the network administrator.

License Key

The 32 Karat software requires a license to collect and analyze data. Without the license, the software will only operate in Demo mode. In Demo mode, only the data supplied with the software in the “data samples” folder can be analyzed and data acquisition is not possible. For instrument control, data acquisition, and analysis, the 32 Karat license key is required. The license key contains the data system serial number and determines other options that can be used with your system. This key is a USB flash drive that must be put in an available USB slot on the controller.

CAUTION

The flash drive license key must remain in place at all times while the software is running.

The license key can be uninstalled on one controller and installed on another, if desired.


P/ACE™ MDQ Plus Software

Launch the P/ACE™ MDQ Plus software by double-clicking the P/ACE™ MDQ Plus software icon on the Windows desktop. The **32 Karat Software Enterprise** window is displayed.

Enterprise Window

Click **View > Hierarchy Pane**. In the left pane, the group **P/ACE MDQ plus** appears. Select this group to display the available instruments in the right pane ([Figure 5.1](#)).

Figure 5.1 32 Karat Software Enterprise Window

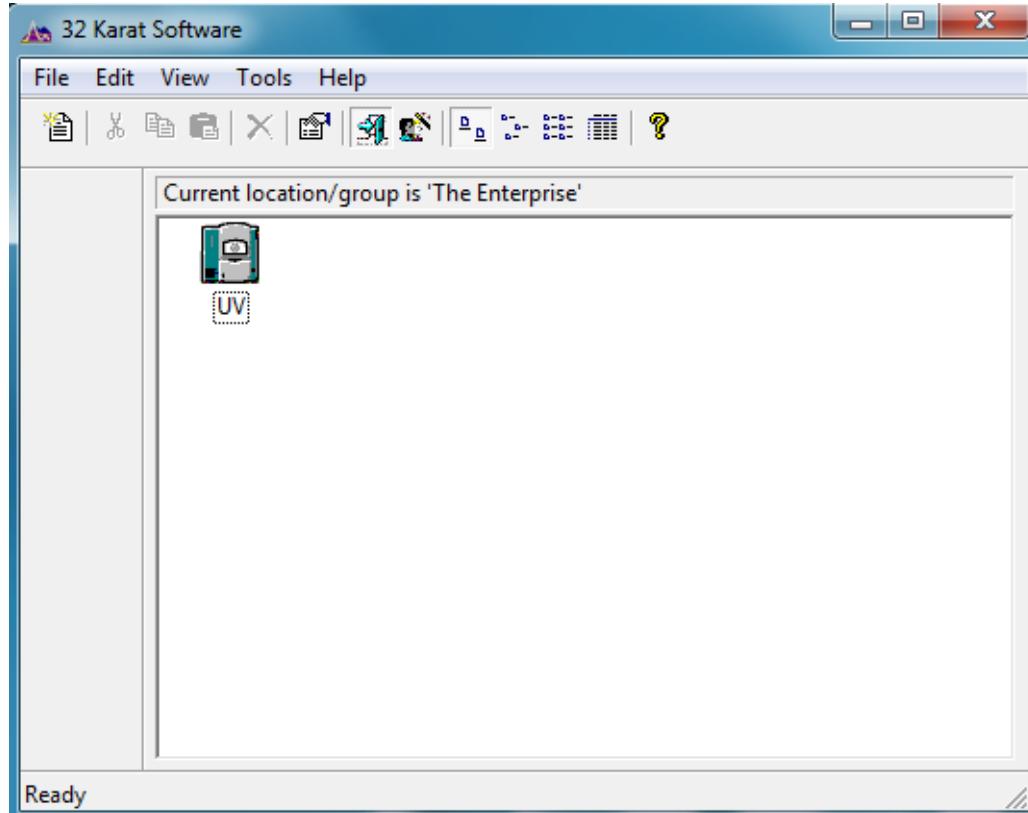Instrument Start Up

While the instrument electronics are off, open the sample and cartridge covers. Verify that a detector and a cartridge are installed. If it is necessary to change the detector, follow the procedure in the maintenance section of the *System Maintenance Guide*. If no cartridge is installed, install the OPCAL cartridge. Check the coolant level for the capillary cooling system by viewing it through the coolant sight glass. If the coolant level is low, add coolant as described in the maintenance section of the *System Maintenance Guide*. Close the cartridge and sample covers and push in the convenience button on the front bezel of the instrument. Allow the instrument up to five minutes to initialize. Shortly after the convenience switch is turned on, the instrument will register a short tone. This tone indicates the circuit boards have activated and started to initialize. Approximately two minutes later, you should hear additional sounds as the filter wheel, transport,

and pressure systems initialize. The instrument is now ready for you to control using the P/ACE™ MDQ Plus software.

System Administration

System administration is a feature used to control access to the software. Access to the P/ACE™ MDQ Plus instrument and the 32 Karat software can be restricted by user, instrument, or project. At installation, system administration is enabled. The default user name is “mdq” and the default password is “plus”. This user name and password are needed for any system changes or access to any instrument. The administrator can grant access to additional users.


To manage user access, click **Tools > Enterprise Login**. When prompted, type the user name, *mdq*, and the password *plus* ([Figure 5.2](#)).

Figure 5.2 System Administration Login Dialog

When system administration mode is enabled, click **Tools > Options**. Select the **Enterprise** tab from the **Options** dialog ([Figure 5.3](#)).

Figure 5.3 Enterprise Tab in the Options Dialog

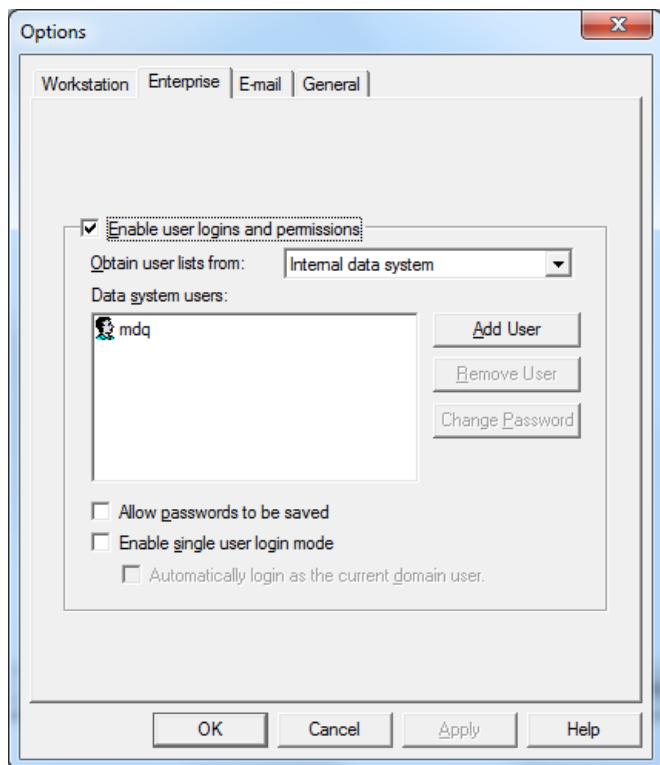
When a user is added as a system user, the access for this user is customized using the System Administration wizard. To access the system administration wizard, click **Tools > System Administration Wizard**.

A project is a set of folders for file types commonly used in P/ACE™ MDQ Plus software. For example, if you are doing analysis with optional LIF detection and another analysis with UV detection, it would be convenient to store your method and report files in folders with the appropriate names. A project is created for each application for the P/ACE™ MDQ Plus system.

When opening an instrument, the user is asked to specify a project. When opened, all folders default to this project. This simplifies locating method and data files. When users have logged on to an instrument, they can select any project they have access to from the **File** menu ([Figure 5.4](#)).

Figure 5.4 File > Project

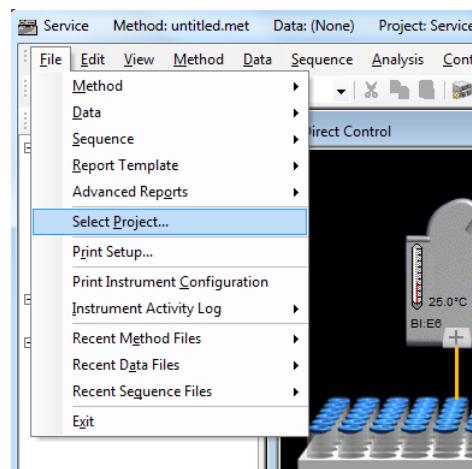Users can lock the **Instrument** window for security if they step away from the instrument. To lock or unlock the window, select lock from the **Window** menu.

Another advantage of using the system administration features is that it stores your customized settings for each instrument. This is convenient when multiple users run different applications on the same system. System administration enables a much greater degree of security on your system.

Configuration

All of the instrument configurations required to run P/ACE™ MDQ Plus system applications are user-defined.

Online versus Offline


The 32 Karat software has two main modes of operation:

- Online mode gives full instrument control and the ability to collect and analyze the data.
- Offline mode allows full data analysis but no instrument control.

If an instrument is running and collecting data, it is possible to open a second window for the same instrument offline. Similarly, if you need to reprocess the data on a remote computer, you should open the instrument offline. The instrument is opened online by double-clicking the instrument icon. To open the instrument offline, right-click the instrument icon and select **Open Offline**. For the purpose of basic training, open the Performance instrument online with the Performance project offline.

When a newly configured instrument is launched for the first time, the **Instrument** wizard appears ([Figure 5.5](#)).

Figure 5.5 Instrument Wizard

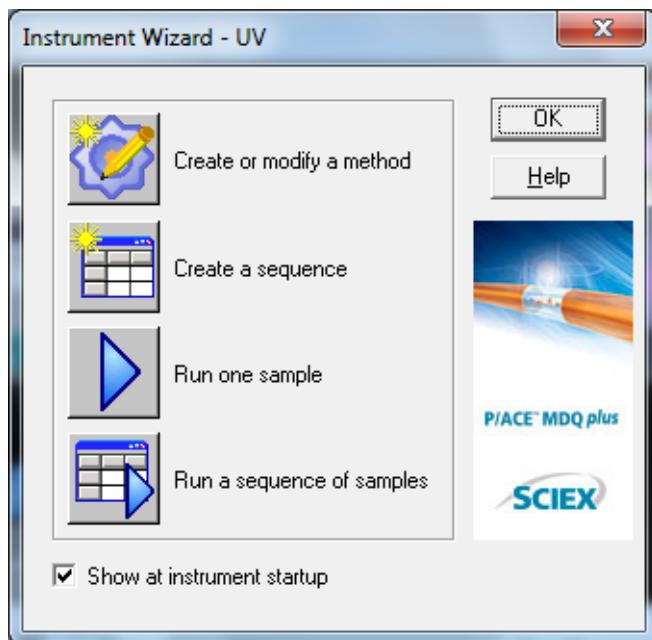
The **Instrument** window normally starts in a small window. Select the square box button in the top-right corner of this window to display the **Instrument** window in full-screen view.

How to Load the Cartridge and Samples


With the instrument and controller initialized, it is good practice to perform some quick checks on the system to make sure performance is correct. In the previous section, the user opened the Performance instrument online. The current project should be “Performance”. If the instrument is open under a different project, it is possible to select the Performance Tests project from the **File** menu on the **Instrument** window ([Figure 5.6](#)).

Figure 5.6 File > Project

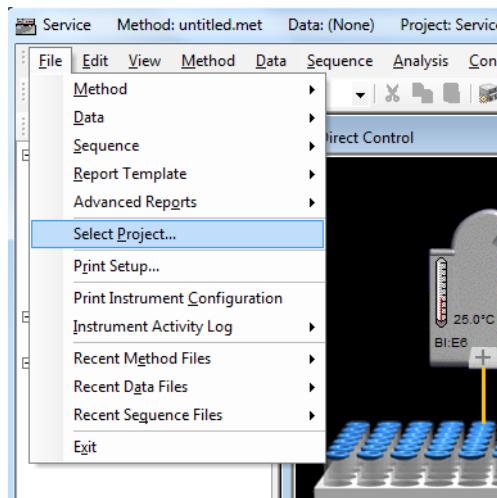
Click **Direct Control** ([Figure 5.7](#)).

Figure 5.7 Control > Direct Control

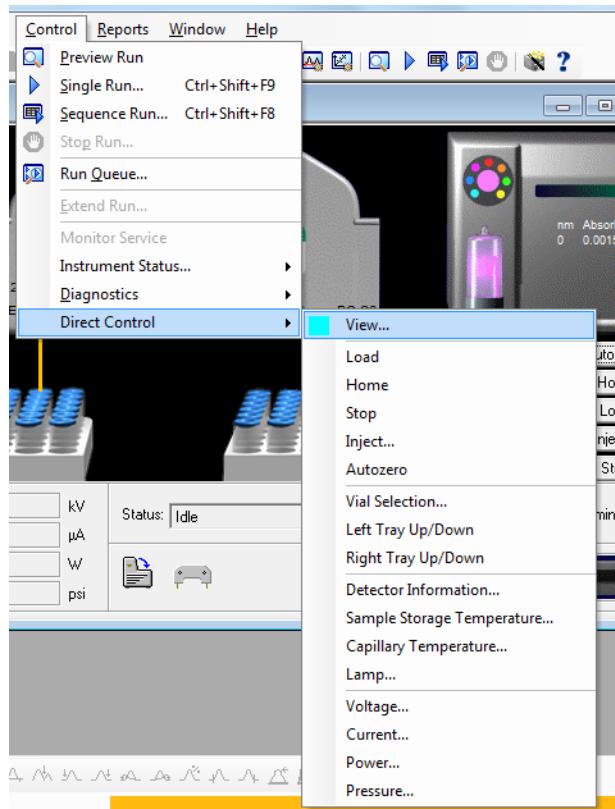The **Direct Control** window allows manual control and reports the status of the instrument functions. Icons that confirm the status of the sample cover and cartridge are located in the lower part of the window ([Figure 5.8](#)). Verify the indicators show the sample cover closed and a cartridge installed. If the sample cover is open, or if no cartridge is installed, the trays cannot move.

Figure 5.8 Direct Control Window

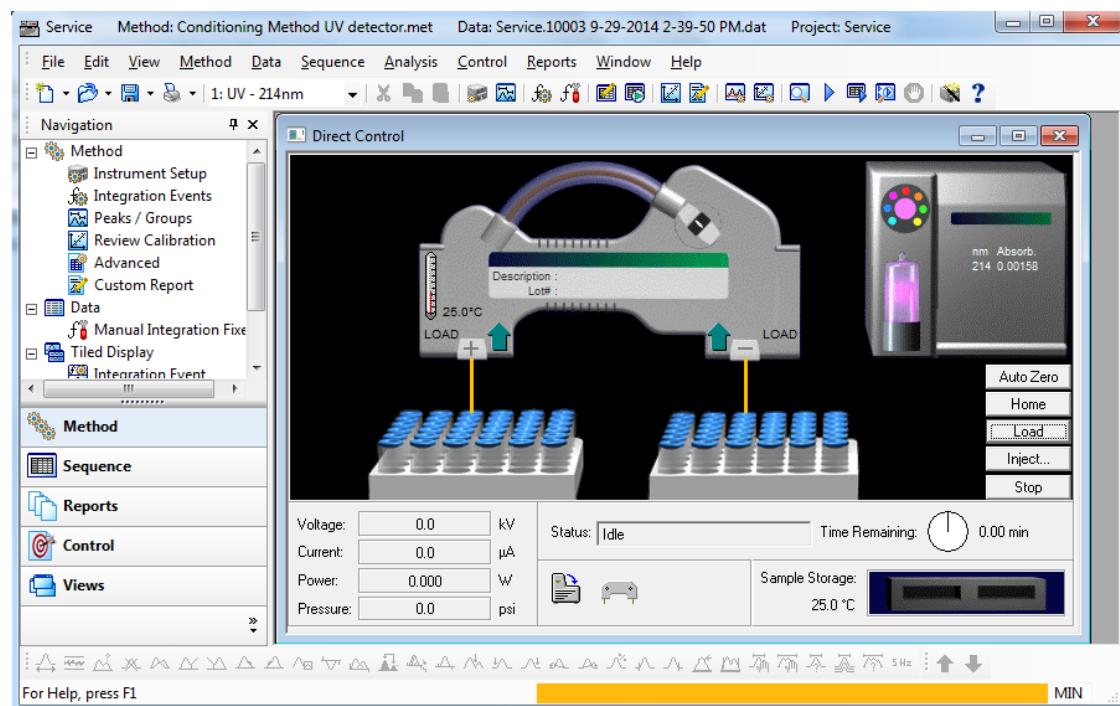Click **Load** to move the trays to the front of the sample compartment. Inspect the interface block by using an angled mirror and penlight. Clean the interface block as required using the procedure in the *System Maintenance Guide*.

Before you run any application, the correct capillary and vials must be installed in the instrument. Install the cartridge and vials required to run the selected method and close the sample cover. The instrument automatically detects the trays and moves to a home position in relation to the interface block. It is now ready to run the method.

How to Run a Method

In the previous sections, the following steps have been taken to prepare the instrument to run a method:

- Start the controller and launch the 32 Karat software.
- Verify the instrument has a detector installed.
- Verify the instrument has a cartridge installed
- Turn on the instrument electronics by engaging the convenience switch.
- Open an instrument online.
- Install a cartridge with the capillary as required by the selected application.
- Install the correct vials and trays for the application.


Once these steps are completed, a method can be run.

-
- 1 Click Single Run or Control > Single Run on the Instrument window (Figure 5.10).

Figure 5.9 Single Run Icon


Figure 5.10 Control > Single Run

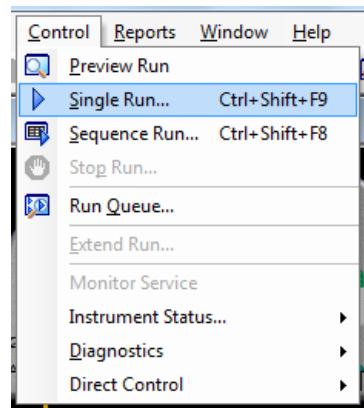
- 2 The Single Run Acquisition dialog appears (Figure 5.11). Click the folder icon on the right of the Method field and select a method. A unique data file name must be entered in the data file field. The data file path will default to the data folder in the currently selected project. Optionally, add additional information for sample description.

Figure 5.11 Single Run Acquisition Dialog

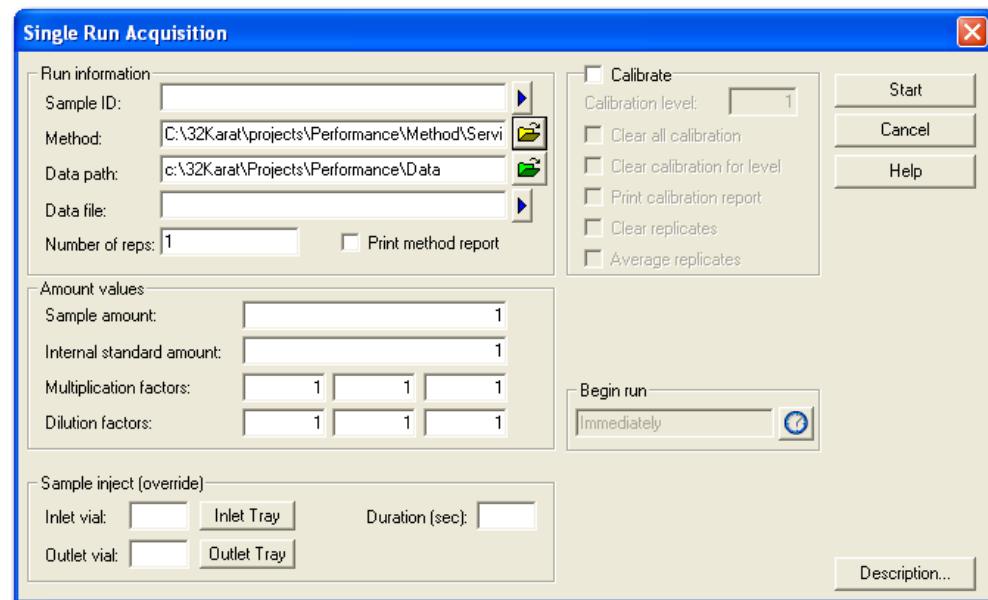
- 3 The instrument runs the selected method. It is not necessary to wait for the method to complete before you start another run. Use the same steps to start a single run. When the run is submitted, it is added to the run queue ([Figure 5.12](#)). Notice the **Submit** button instead of **Start**.

Figure 5.12 Single Run Acquisition Dialog—Running a Method

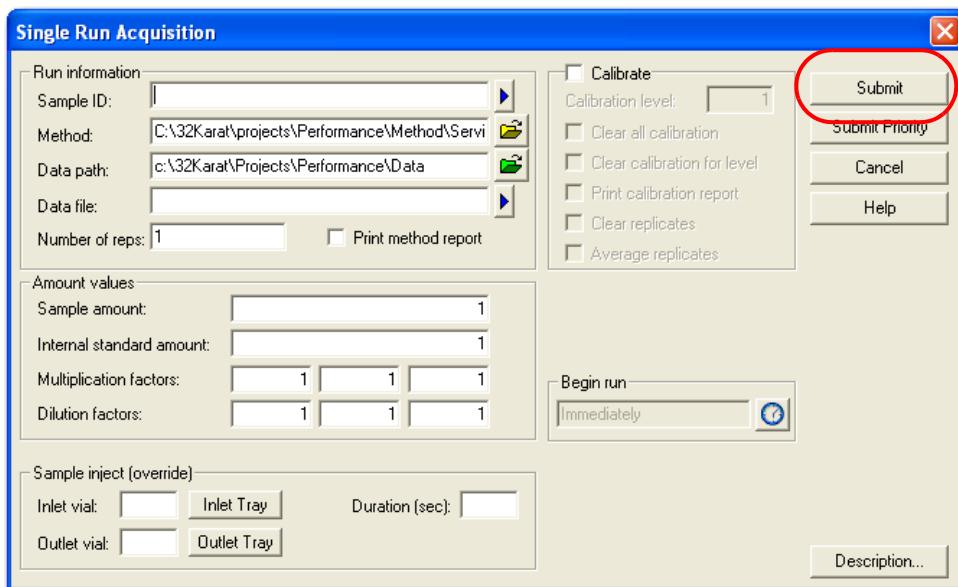- 4 To view the run queue, click **Control > Run Queue** ([Figure 5.13](#)).

Figure 5.13 Control > Run Queue

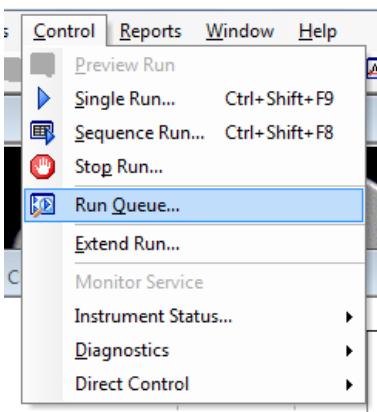
Analysis

After the data is created using a method, it must be analyzed.

- 1 Click **File > Data > Open**.
- 2 Open the data file collected with the method that was run as described in the previous section.
- 3 Click the **Window** menu at the top of the **Instrument** window.
- 4 At the bottom of the **Window** menu, the available data channel views appear. Click **UV - 214nm**. An electropherogram with two main peaks appears.

The menu on the electropherogram should display “1:UV - 214nm”. This menu indicates which channel of data that the selected analysis and window settings will change. Select from the menu to display other channels when present ([Figure 5.14](#)).

Figure 5.14 Analysis Channel Menu

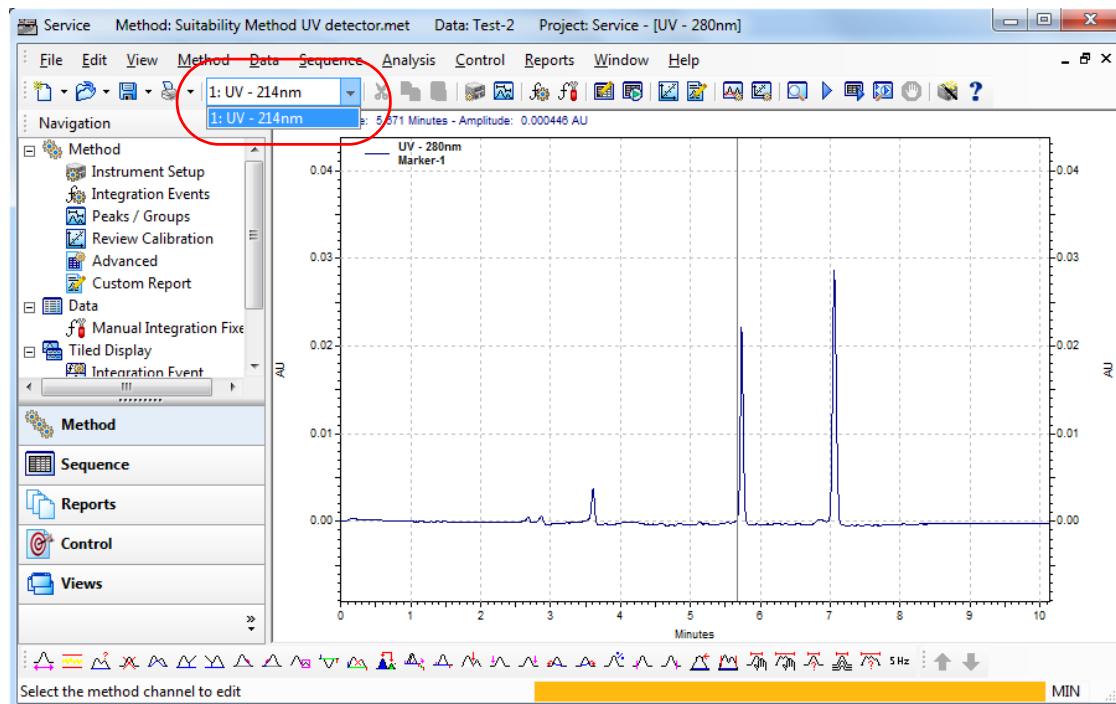- 5 Click **Channel 1 (1:UV - 214nm)** for analysis and return to the electropherogram window.

NOTE An electropherogram is a plot with migration time on the X axis and absorbance data plotted on the Y axis. These plots generally have long sections of relatively flat detector response known as baseline. If a compound is detected, a large deflection is shown from the baseline, known as peaks. You must accurately identify what parts of the electropherogram are baseline segments and what parts of the electropherogram are peaks. Integration is used to automate this process.

-
- 6 Click the **Integration Events** icon (Figure 5.15) to display the **Integration Events** table (Figure 5.16).

Figure 5.15 Integration Events Icon


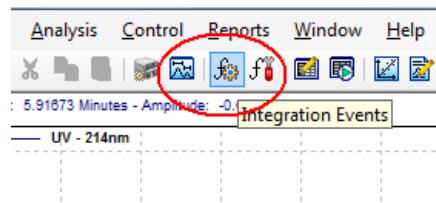
Figure 5.16 Integration Events Table

A screenshot of a table window titled 'Integration Events -- UV - 214nm'. The table has columns for '#', 'Event', 'Start Time', 'Stop Time', and 'Value'. There are four rows, each with a checked checkbox in the first column. Row 1: Event 'Threshold', Start Time 0.000, Stop Time 0.000, Value 2000. Row 2: Event 'Width', Start Time 0.000, Stop Time 0.000, Value 0.2. Row 3: Event 'Shoulder Sensitivity', Start Time 0.000, Stop Time 0.000, Value 10000. Row 4: An empty row with a checked checkbox.

-
- 7 As windows are opened, other windows needed are covered. To arrange the windows, click **Window > Tile Horizontally**. The open windows fill the screen as equal-size tiles. Close all windows except for the **Integration Events** table and the **UV - 214nm** window. Click **Window > Tile Horizontally** so that these two windows fill the screen. The title bar of the **Integration Events** table indicates this screen is for 214 nm. If another channel is indicated, change the analysis channel as described in previously. By arranging the **Integration Events** table and the **UV - 214nm** window, users can observe how changes in the table affect how peaks are detected in the electropherogram (Figure 5.17).

Figure 5.17 UV - 214nm Window and Integration Events Table, Tiled

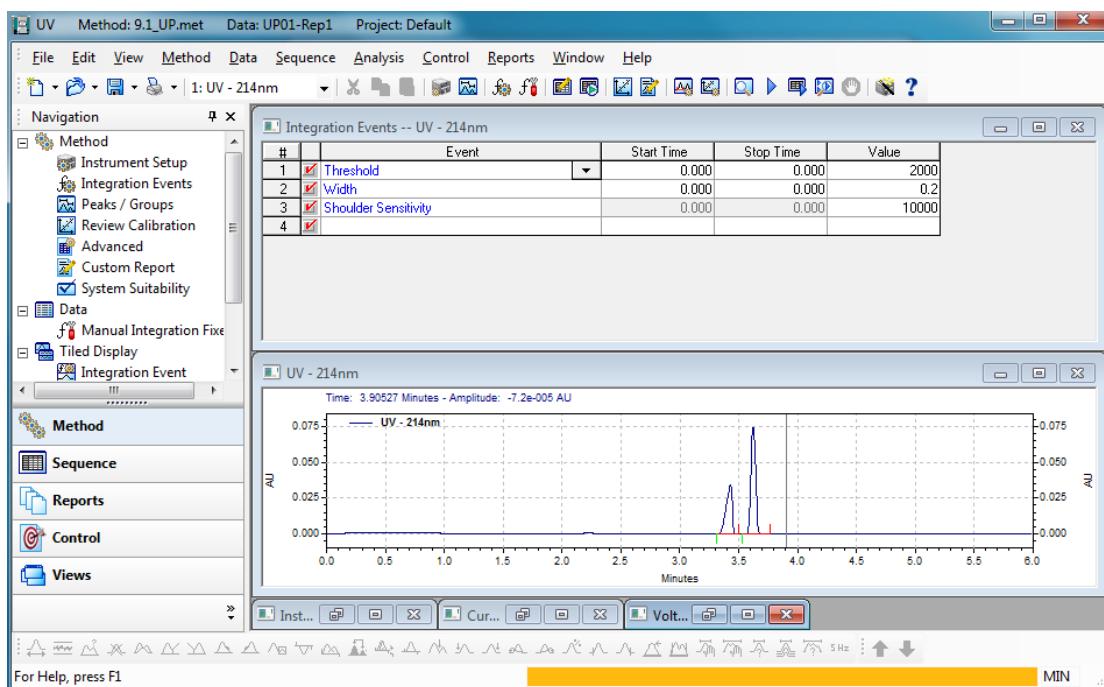8 There are four parameters that require defaults for integration to work:


- **Width**—Adjusts how sensitive peak-detection is for changes in the baseline. Should be set to the width of the smallest peak in the electropherogram. The default value is 0.2.
- **Threshold**—Adjusts the distance from which a point is considered a cluster point. The software considers points that are removed from the baseline by a large enough distance as cluster points. You can picture the threshold as an imaginary line that runs parallel to the baseline above and below it. If a sufficient number of cluster points are grouped together, the software calls this group of points a peak. A curve is drawn through these points that extends to the baseline. The points where the curve meets the baseline are the start and stop points.
- **Shoulder Sensitivity**—Measures the curvature of the upslope and downslope of the peak.
- **Minimum Cluster Distance**—How many data points must be grouped together before the software identifies the group of cluster points as a peak.

Optimization of the integration parameters is described in the *Methods Development Guide*.

Area Calculation

In addition to the calculation of baseline, peak start, peak apex, and peak stop, the integration algorithm is responsible for calculating the area under the peak. The area calculation is a simple numeric calculation. The area under the peak is calculated for each data point. All of the areas for the data points of a peak are added together to calculate the area of the peak. For a given method, if the same sample is injected several times, the peaks should occur at the same times and corresponding peaks should have the same areas. The time of the peak relates to the identification of the compound while the area of the peak relates to the quantity of the compound. These two parameters are not completely independent. For example, if a longer pressure injection is used, a larger amount of sample is introduced. This sample is injected farther into the capillary and as a result, when voltage is applied, it reaches the detector earlier. The peak area will be larger because there is more sample present. [Figure 5.18](#) is an example of how injection volume can change migration time and peak area.

Figure 5.18 Injection Volume vs. Migration Time and Peak Area

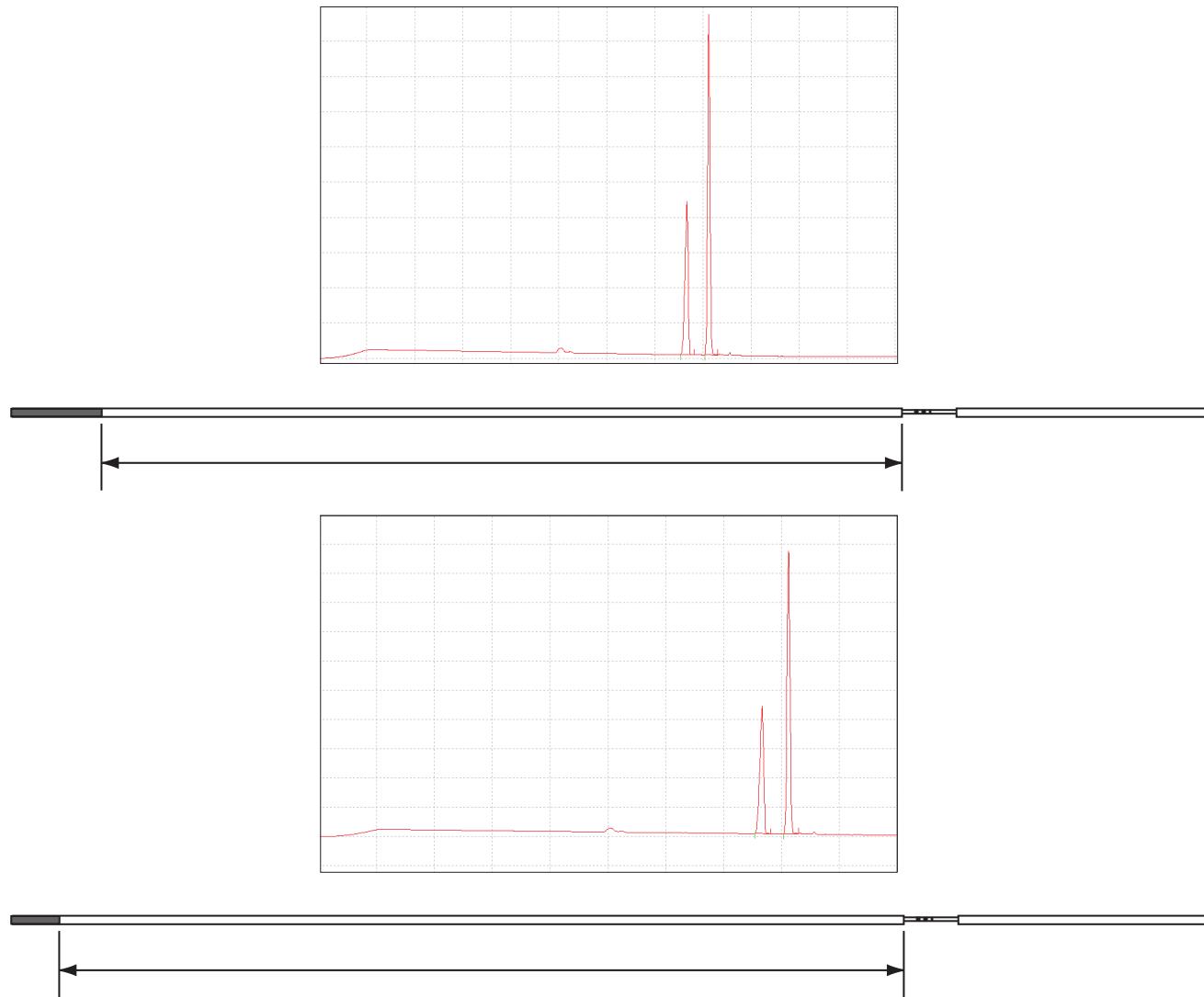
Capillary condition or temperature also affects these parameters. A cooler temperature causes the viscosity of the solutions to change. This in turn causes a change in injection volume and the parameters change as previously described. In these two examples, the performance of the instrument is responsible for the change. If the capillary condition is not good, the migration time can change due to a condition known as electroosmotic flow. This condition is the result of ionization of the capillary wall. This ionization causes a bulk fluid flow through the capillary from the positive to negative electrode. This flow adds or subtracts to the velocity of the compound's inherent mobility.

The compound reaches the detector at a different time because it moves at a different velocity. The change in velocity causes the compound to stay in front of the detector for a different amount of time. This causes the detector response to be different. A slow moving peak is reported as a larger peak than a similar sized fast-moving peak. This condition is not caused by instrument performance, but is minimized in an optimal method. To further adjust for this condition, the peak area reported with the peak velocity can be corrected. The resultant area value is called corrected area. The 32 Karat software can calculate this amount automatically.

Annotations

Select the box near the upper-right corner of the **UV - 214nm** window to maximize it. If the data is analyzed, baseline and peak marks are present on the electropherogram. The baseline segments are shown in red across the bottom of the peaks. Short vertical tick marks indicate the peak start and stop locations. The peak start is drawn in red below the baseline segment and the peak stop is drawn in blue above the baseline segment. Right-click the display and select **Annotations**. The **Trace Annotations Properties** dialog will display. At the bottom of the dialog are options for general graph properties.

Figure 5.19 Trace Annotation Properties Dialog

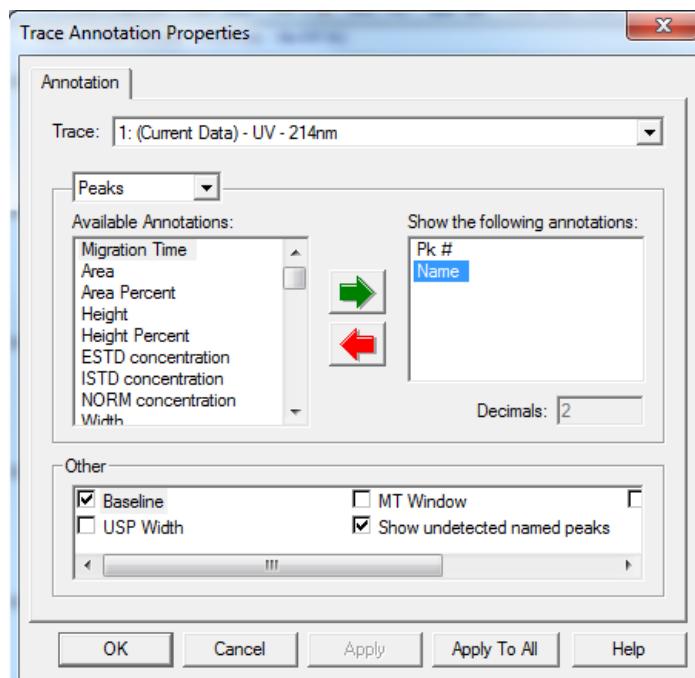By default, the **Baseline** and **Show undetected named peaks** options are selected. The **Available Annotations** option contains parameters that can be displayed for each peak. Click **Peak #** and **Name** using the green arrow to move them to the right pane (Figure 5.19) and then click **Apply**. As a result the peak numbers are displayed above the peaks. Peak names are displayed on the electropherogram because they are programmed into the method. Depending on the quality of the separation, the peak names may not be associated with the correct peaks at this point in time.

Peak Identity


Add information to the **Peak** table to name a peak. This information is already entered in the **Peak** table for the default methods supplied with the P/ACE™ MDQ Plus system. In other supported applications with unknown samples, it may be necessary to add new peak information to the **Peak** table. For details on adding peaks to the **Peak** table, refer to the *Methods Development Guide*.

Run Multiple Methods with a Sequence

Once the method and analysis are optimized, the method can be run repeatedly on different samples or multiple times on a single sample to calculate the performance. To automate this process, program a sequence. The easiest way to program a sequence is by using the **Sequence** wizard. The **Sequence** wizard runs the same method the specified number of times. There are many additional selections for run customization using the **Sequence** wizard. For details on using the **Sequence** wizard, refer to the *Methods Development Guide*.

- 1 Click **Sequence Process**.

Figure 5.20 Process Sequence Dialog

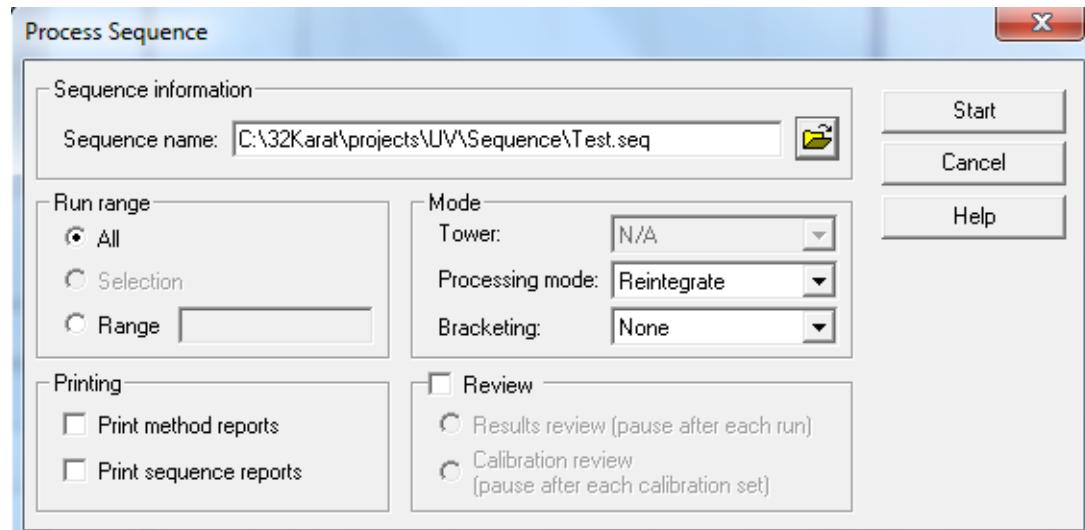- 2 Click the folder icon.

-
- 3 Click the sequence.
 - 4 In **Run range**, specify whether to process the complete sequence (**All**), or only a number of runs (select **Range**).

To specify a number of runs, type the sequence line numbers in the **Range** text box. Separate ranges by using dashes and commas. For example, to run lines 2, 3, 4, and 6, type the **Range** as **2-4, 6**.
 - 5 Click **Start** to process the sequence.
-

Observe that as the sequence is processed the status of each line is updated. When the run is completed, the last data analyzed will display in the **Instrument** window. The analysis of each data file is shown only for a moment until the next line of the sequence is processed. To use a report, the run type must be changed in the sequence. Click **Sequence > Edit** to open the sequence. Click and drag over the **Run Type** cells to select the cells. Right-click and click **Set Run Type**. Click **Summary Run** as the **Run Type**.


When **Summary Run** is selected, the text in the Run Type cells changes. The first cell is “Summary Begin”, the last cell is “Summary End”, and the other cells are “Summary Types” ([Figure 5.21](#)). When the sequence is processed a summary report generates.

Figure 5.21 Sequence Window

Run #	Status	Run Type	Level	Conc Override
1		Summary Begin	0	n/a
2		Summary Run	0	n/a
3		Summary Run	0	n/a
4		Summary Run	0	n/a
5		Summary End	0	n/a
6				


To specify which summary is reported, double-click the **Summary Begin** cell in the **Run Type** column. Observe that the **Begin Summary** check box is selected ([Figure 5.22](#)). Click the text next to the check box. In the right side of the window a report template can be specified. Select a report template file and click **OK**.

Figure 5.22 Sample Run Type(s) Dialog

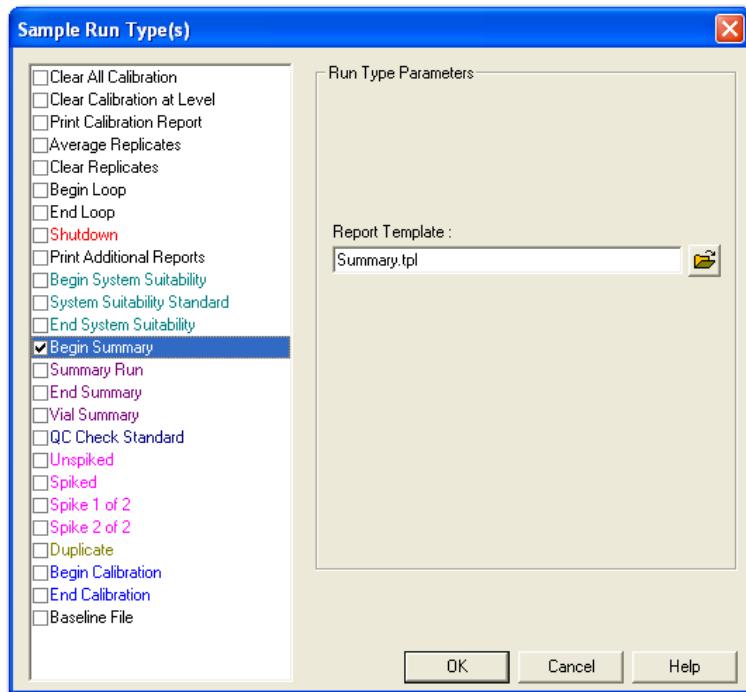
Reprocess the sequence. Select **Reports**. All available sequence reports display. The example sequence only has one report specified, open the displayed report ([Figure 5.23](#)).

Figure 5.23 Reports > View > Sequence Report

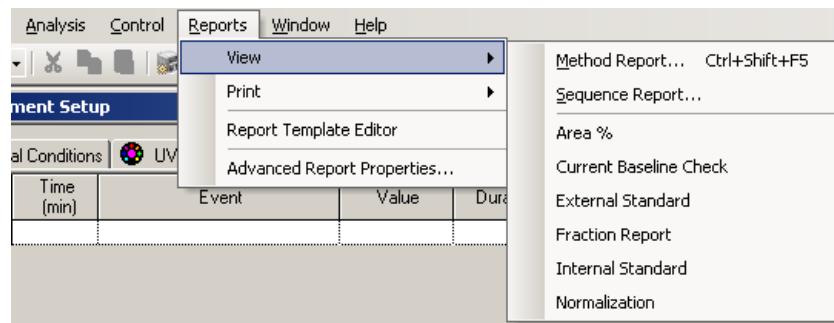
When a sequence is created for acquisition, the **Sequence** wizard provides more options. Use the **Sequence** wizard to create a sequence for acquisition with the same method, instead of selecting data files from a directory, give a data file name. When the sequence is run, the collected data is saved under this name. To prevent over-writing data, the file can be appended with one or more fields. The available fields are shown by selecting the arrow to the right of the filename field ([Figure 5.24](#)). Click **Next**.

Figure 5.24 Sequence Wizard—Data Path Menu

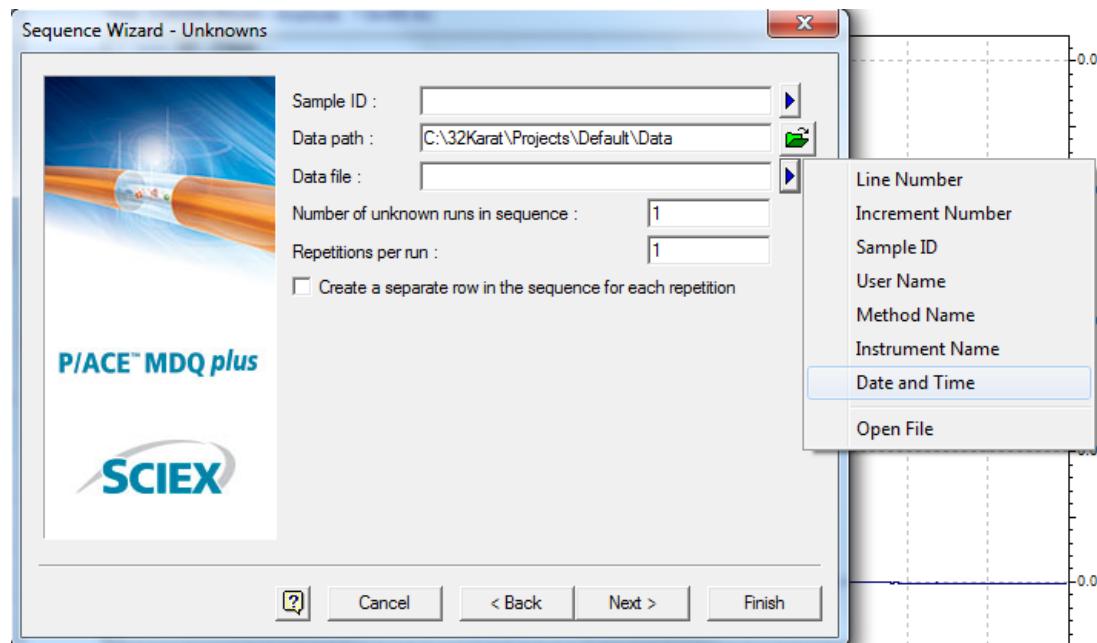In the **Vials** window, the positions of the sample vials are specified. Vials can be incremented by row or by column (Figure 5.25). Click **Next >**.

Figure 5.25 Sequence Wizard—Vials

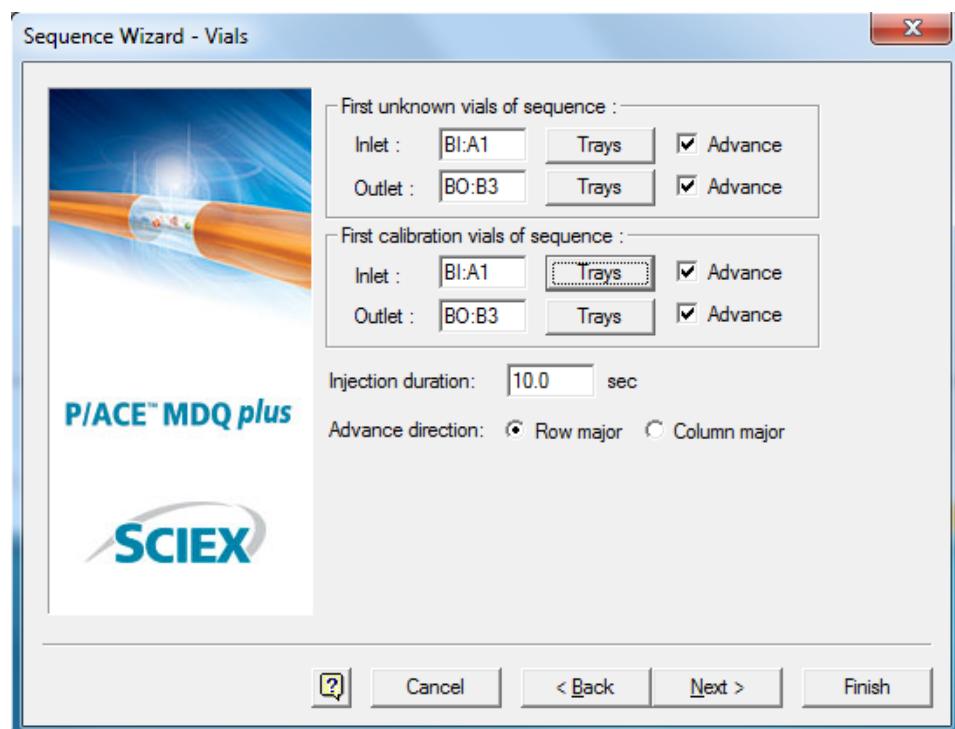Viewing Reports

Standard Reports

Each project contains certain default method reports. After each data analysis, parameters for certain peaks will be reported. For the area% report, all detected peaks are included. For ESTD, ISTD, and normalization reports, only detected peaks included in the peak table are reported. These reports cannot be edited.

Method Custom Reports

Method custom reports can be modified to display custom information and saved as a new report. When you include a method custom report with a method, a method custom report is generated every time data is analyzed.

-
- 1 Click **Method Custom** or **File > Report Template Open**.

 - 2 Select the report template, for example the external standard report template, **External Standard.Srp**.

 - 3 Modify the report by right-clicking on any section of the report.

 - 4 Click **File > Report Template > Save as** and save the report as *CRPtrain.srp*.
-

Sequence Custom Report

On the surface, sequence custom reports look very similar to the custom method reports. However, the available parameters are more particular to the sequence than to the data. These reports only generate after processing a sequence, if the report template is specified as part of the sequence.

Advanced Reports

Detailed information on advanced reports available from the *32 Karat Software Online Help* that is delivered with the P/ACE™ MDQ Plus system.

-
- 1 Click **File > Advanced Reports > New**.

 - 2 Right-click cell **A1** of the report window and then click **Table Wizard**.
-

-
- 3 Click **Sequence Summary Table** and then click **Next**.
A list of peak parameters appears.
- 4 To obtain a sequence summary of Area and Migration Time, select those parameters and click **Add** to add them to the right panel of the window.
For this example, click **Area** and **Migration Time**.
- 5 The **Index** is used to select which trace is reported on. When Index is selected, a list of traces appears.
- 6 The **Precision** field, located at the bottom of the window, determines the number of digits after the decimal point. Click **Next**.
- 7 The current window contains options for the types of peak in the report. Click option 1 name and option 2 name to report all named peaks.
- 8 On the current window specify the data file parameters by clicking the data file name and selecting the right arrow. Click **Next**.
- 9 The current window configures the data display on reports (by row or by column). Click **Next**.
- 10 This window asks if you want the report to display data by row or by column. Click **Next**.
- 11 The current window contains the option for displaying statistics. Click **Yes** and **Finish**.
- 12 Click **File > Advanced reports > Save as** and save the report as a template using the file name *RSD.tpl*.
This file is a report shell that is ready to summarize and report a sequence when it is run.

Suitability Reports

A suitability report is a special type of sequence report. To generate a suitability report, create a sequence with Begin Suitability and End Suitability run types. From the **Begin Suitability Run Type** window, specify the template. To change the suitability parameters included in the report, click **Method** from the **Instrument** menu.

NOTE System Suitability is only available if System Suitability is selected in the instrument configuration. All default P/ACE™ MDQ Plus instrument configurations have System Suitability enabled.

Viewing and Printing Reports

To view or print a report, click **Reports** on the **Instrument** window. If the **Print Reports** option is selected on the **Run** window, method or sequence reports are automatically printed.

Maintenance

Refer to the *System Maintenance Guide* for the following procedures.

- Install the UV or PDA Detector
- Install Wavelength Filters for the UV Detector
- Calibrate the PDA Detector
- Install the LIF Detector
- Calibrate the LIF Detector
- Install the Capillary Cartridge
- Install the Universal Vials and Caps
- Clean the Ejecting Levers, Electrodes, and Interface Block
- Refill the Capillary Coolant
- Replace the Deuterium Lamp
- Procedures for Instrument Care
- Procedures for LIF Detector Care and Maintenance

Specifications

P/ACE™ MDQ Plus System

Table A.1 System Specifications

Item	Description
Dimensions (H x W x D)	99.1 cm (cover open), 73.7 cm (cover closed) × 63.5 cm × 72.4 cm (39 inches (cover open), 29 inches (cover closed) × 25 inches × 28.5 inches)
Weight	85.3 kg (188 lbs)
Electrical	<p>Power requirement: 100 VAC to 240 VAC, 5.0 A, 50 Hz or 60 Hz</p> <p>Power consumption: supply voltage must not exceed 10% of nominal</p> <p>Fuses (depending on supply voltage in use):</p> <ul style="list-style-type: none"> • 8.0 A slow blow; 1/4 inch (2 ea.): 100 VAC to 120 VAC • 6.3 A time delay; 20 mm (2 ea.): 200 VAC to 240 VAC <p>Installation (overvoltage) category: Category II</p>
Working Environment	<p>Altitude: up to 2000 m (6562 feet)</p> <p>Humidity: <80% (non-condensing) at 15 °C to 30 °C <60% (non-condensing) at 30 °C to 40 °C</p> <p>Temperature: 15 °C to 40 °C (15 °C to 30° C recommended)</p>
Maximum Heat Dissipation	400 W (1024 BTU/hour)
Pollution Degree	2
I/O	<p>TTL: 2</p> <p>Contact closures: 2</p>

Validated Controller Configuration

The system includes a computer and a monitor (also referred to as a “controller”). The software has been validated on a controller with the following specifications:

- Lenovo M720s workstation
- 8th Generation Intel Core i5-8600 3.1 GHz processor
- 8 GB RAM
- 22 inch wide-screen monitor with True Color and 1680 x 1050 resolution
- Windows 10 Enterprise LTSB 2016 (Windows 10 IoT) with Cybersecurity
- Operating system language set to English (United States)
- 500 GB hard drive
- DVD-RW drive
- 2 serial ports
- 2 Ethernet ports
- 8 USB ports

NOTE SCIEX fully validates and supports the controllers supplied with the system. Only limited support is available for customer-supplied computers.

NOTE Specifications are subject to change without notice.

Sample Temperature Control

Table A.2 Sample Temperature Control Specifications

Specification Type	Description
Temperature Range	20 °C below ambient to 60 °C (140 °F) Minimum setting: 4 °C (39 °F)
Temperature Stability	±1 °C at 25 °C (77 °F) ±3 °C at 4 °C (37 °F) and 60 °C (140 °F)
Temperature Accuracy	±2 °C within a range of ±15 °C from ambient temperature ±3 °C outside a range of ±15 °C from ambient temperature

Capillary Temperature Control

Table A.3 Capillary Temperature Control Specifications

Specification Type	Description
Temperature Range	10 °C below ambient to 60 °C (140 °F) Minimum setting: 15 °C (59 °F)

Table A.3 Capillary Temperature Control Specifications (*Continued*)

Specification Type	Description
Temperature Stability	±1 °C at 25 °C
Temperature Accuracy	±1 °C within a range of ±1 °C from ambient ±2 °C outside a range of ±5 °C from ambient

Pressure and Vacuum System

Table A.4 Pressure and Vacuum System Specifications

Specification Type	Description
Pressure Range	Injection: 0.1 psig to 25 psig (pressure) or 0.1 psig to 5.0 psig (vacuum) Rinse: 0.1 psig to 100 psig (pressure) or 0.1 psig to 5.0 psig (vacuum)
Pressure Stability	±0.3 psi at 25 psi ±1.0 psi at 100 psi
Pressure Direction	Applied at inlet or outlet for all pressure functions, rinses, and injections User-settable in software

Detector Specifications

UV Detector

Table A.5 UV Detector Specifications

Specification Type	Description
Wavelength Range	190 nm to 600 nm
Wavelength Accuracy	±2 nm
UV Source	30 W pre-aligned deuterium lamp
Filter Selection	200 nm, 214 nm, 254 nm, and 280 nm (standard) with 3 open positions for additional filters Filter dimensions must be: <ul style="list-style-type: none">• Diameter: ½ inch (127 mm)• Thickness: 0.20 inch (5 mm)
Analog Output	Output 1 is data Full scale output is 1.0 AU/V. Multipliers of 1.0, 0.5, 0.2, 0.05, 0.02, and 0.01 to provide lower AU/V values can be set in the software. Output 2 not used Output 3 depends on what is programmed: <ul style="list-style-type: none">• Current signal when voltage is programmed• Voltage signal when current or power is programmed

Laser Induced Fluorescence (LIF) Detector (Optional)

Table A.6 LIF Detector Specifications

Specification Type	Description
Wavelength Range (for optics)	Excitation: 300 nm to 700 nm Emission: 350 nm to 750 nm
Solid State Laser	Laser output delivered to the capillary: 2.5 ± 0.5 mW Laser wavelength 488 nm nominal
Sensitivity	Minimum signal/peak-peak noise ratio of 10,000:1 for 50 nM sodium fluorescein in a 75 μm i.d. capillary
Relative Fluorescence Units Range	0 RFU to 1000 RFU
Filters (optional)	For 488 nm laser: 488 notch filter and 520 nm band-pass filter For user-supplied lasers, two filters are required: a filter to block stray laser light and an emission filter to select the wavelength of the emitted light. Filter dimensions must be: <ul style="list-style-type: none"> • Outer diameter: 0.500 inch (+0.000 inch, -0.010 inch); 12.7 mm (+0.000 mm, -0.25 mm) • Thickness: ≤ 0.350 inches (0.889 mm) • For multiple filters used in a single channel, total thickness: ≤ 0.350 inches (0.889 mm)
Dynamic Range	$> 10^4$
Baseline Noise	< 0.005 RFU peak to peak
Baseline Drift	< 0.2 RFU/hour
Analog Outputs	Output 1 is Data Channel 1 Output 2 is Data Channel 2 Full scale output is 1.0 AU/V. Multipliers of 1.0, 0.5, 0.2, 0.05, 0.02, and 0.01 to provide lower AU/V values can be set in the software. Output 3 depends on what is programmed: <ul style="list-style-type: none"> • Current signal when voltage is programmed • Voltage signal when current or power is programmed

Photo Diode Array (PDA) Detector (Optional)

Table A.7 PDA Detector Specifications

Specification Type	Description
Wavelength Range	190 nm to 600 nm
Wavelength Accuracy	2 nm
UV Source Lifetime	1000 hours
UV Source	30 W pre-aligned deuterium lamp
Scan Collection Frequency	0.5 Hz to 32 Hz
Detector	256 element diode array
Bandwidth	6 nm minimum (absorbance averaging)

Table A.7 PDA Detector Specifications (*Continued*)

Specification Type	Description
Analog Output	<p>Output 1 is Data Channel 1</p> <p>Output 2 is Data Channel 2</p> <p>Full scale output is 1.0 AU/V. Multipliers of 1.0, 0.5, 0.2, 0.05, 0.02, and 0.01 to provide lower AU/V values can be set in the software.</p> <p>Output 3 depends on what is programmed:</p> <ul style="list-style-type: none">• Current signal when voltage is programmed• Voltage signal when current or power is programmed

Customer Training

- In North America: NA.CustomerTraining@sciex.com
- In Europe: Europe.CustomerTraining@sciex.com
- Outside the EU and North America, visit sciex.com/education for contact information.

Online Learning Center

- [SCIEXUniversity](#)

SCIEX Support

SCIEX and its representatives maintain a staff of fully-trained service and technical specialists located throughout the world. They can answer questions about the system or any technical issues that might arise. For more information, visit the SCIEX website at sciex.com or contact us in one of the following ways.

- sciex.com/contact-us
- sciex.com/request-support

CyberSecurity

For the latest guidance on cybersecurity for SCIEX products, visit sciex.com/productsecurity.

Documentation

This version of the document supersedes all previous versions of this document.

To view this document electronically, Adobe Acrobat Reader is required. To download the latest version, go to <https://get.adobe.com/reader>.

To find software product documentation, refer to the release notes or software installation guide that comes with the software. Documentation for the hardware products can be found on the Customer Reference DVD that comes with the system or component.

For the latest versions of the documentation, visit the SCIEX website at sciex.com.

NOTE To request a free, printed version of this document, contact sciex.com/contact-us.